



An optimized toolkit for high molecular weight DNA extraction and ultra-long-read nanopore sequencing using glass beads and hexamminecobalt(III) chloride

Inswasti Cahyani, John Tyson, Nadine Holmes, et al.

Genome Res. published online March 24, 2025

Access the most recent version at doi:[10.1101/gr.279943.124](https://doi.org/10.1101/gr.279943.124)

P<P Published online March 24, 2025 in advance of the print journal.

Accepted Manuscript Peer-reviewed and accepted for publication but not copyedited or typeset; accepted manuscript is likely to differ from the final, published version.

Open Access Freely available online through the *Genome Research* Open Access option.

Creative Commons License This manuscript is Open Access. This article, published in *Genome Research*, is available under a Creative Commons License (Attribution-NonCommercial 4.0 International license), as described at <http://creativecommons.org/licenses/by-nc/4.0/>.

Email Alerting Service Receive free email alerts when new articles cite this article - sign up in the box at the top right corner of the article or [click here](#).

Comprehensive immune receptor profiling.
Discover the **DriverMap™ AIR Assay** difference.

LEARN
MORE



To subscribe to *Genome Research* go to:

<https://genome.cshlp.org/subscriptions>

Published by Cold Spring Harbor Laboratory Press

An optimised toolkit for high molecular weight DNA extraction and ultra-long-read nanopore sequencing using glass beads and hexamminecobalt(III) chloride

Inswasti Cahyani¹, John Tyson², Nadine Holmes¹, Josh Quick³, Chris Moore¹, Nick Loman³,
Matthew Loose¹

¹School of Life Sciences, University of Nottingham NG7 2UH, United Kingdom

²UBC, Dept of Pathology and Laboratory Medicine, Vancouver, Canada

³School of Biosciences, University of Birmingham B15 2TT, United Kingdom

Corresponding author: matt.loose@nottingham.ac.uk

Abstract

Since the advent of long read sequencing, achieving longer read lengths has been a key goal for many users. Ultra-long read sets ($N50 \geq 100$ kb) produced from Oxford Nanopore sequencers have improved genome assemblies in recent years. However, despite progress in extraction protocols and library preparation methods, ultra-long sequencing remains challenging for many sample types. Here we compare various methods and introduce the FindingNemo protocol that: (1) optimises ultra-high molecular weight (UHMW) DNA extraction and library clean-up by using glass beads and Hexamminecobalt(III) chloride (CoHex), (2) can deliver high ultra-long sequencing yield of >20 Gb of reads from a single MinION flow cell or >100 Gb from PromethION devices (R9.4 to R10.4 pore variants), and (3) is scalable to using fewer input cells or lower DNA amounts, with extraction to sequencing possible in a single working day. By comparison, we demonstrate that this protocol surpasses previous methods by enabling precise determination of input DNA quantity and quality through cell counting, sample dilution, and homogenization techniques.

Introduction

The first ultra-long (UL) sequencing protocols on nanopore sequencing instruments resulted in highly contiguous human genome assemblies (Jain et al. 2018a; Quick 2018). Later methods improved on these assemblies through increased yields (Logsdon et al. 2021; Miga et al. 2020). Oxford Nanopore Technologies (ONT) subsequently launched UL sequencing kits, starting with ULK001 on R9.4 flow cells and more recently ULK114 on R10.4.1 flow cells. The ULK001 protocol utilised extraction and clean-up steps using silicon coated discs from Circulomics Inc. (PacBio

2023), which was subsequently replaced with the ONT star-shaped matrix for library precipitation in the ULK114 kit. These UL kits introduced a range of protocol improvements both in library preparation and sequencing performance to maximise read length and generated reads with an N50 > 100kb from both the MinION and PromethION flow cells.

In our experience, we have found that the most critical step in obtaining high throughput, high occupancy UL sequencing libraries is the DNA fragmentation step, followed by the final clean-up step prior to loading the library. The UL protocol uses a transposase enzyme complex (transposome) to fragment DNA. During fragmentation, the DNA should be homogeneous and accessible to this fast-acting transposome. Two methods are effective in reducing the number of transposase cuts to increase read length: either using a high concentration of high molecular weight (HMW) DNA as shown by Jain et al in sequencing the human genome (Jain et al. 2018a) or diluting the transposase relative to DNA (Jain et al. 2018b). UL libraries are difficult to clean up using conventional magnetic solid phase reversible immobilisation (SPRI) approaches as the DNA can be hard to elute and so sheared during this; SPRI beads are mostly used for PCR clean-up where the fragment length is much shorter (Rodrigue et al. 2010; Stortchevoi et al. 2020). In addition, the library is a complex of motor protein and DNA, and so alcohol-based precipitation and washes are not compatible. In our search for suitable alternative approaches, we found that hexamine cobalt (III) (CoHex) cations stabilise and condense DNA and reasoned they may enable alcohol-free precipitation of nanopore sequencing libraries (Allers and Lichten 2000; Kankia et al. 2001).

Early versions of the ONT ultra-long (UL) protocol required extraction and library clean-up starting from 6 million human cells, or an equivalent input amount of DNA (PacBio 2023, Methods). The complete extraction and library preparation protocol takes between two to three days in the laboratory. However, different sample types often require tailored extractions to provide optimal DNA for sequencing, so we sought to use a wide range of sample extraction protocols to obtain UL reads on ONT platforms. We also explored reducing both sample input requirements and the total time for library preparation for sequencing scalability.

Here we describe the FindingNemo protocol (named after the characteristic orange colour of the CoHex buffer) for the generation of high occupancy ultra-long reads on nanopore platforms. This protocol can generate equivalent or more throughput to disc-based methods and may have

additional advantages in tissues and non-human cell material. The protocol can also be tuned to enable extraction from as few as one million human cell equivalents or 5 µg of human ultra-high molecular weight (UHMW) DNA as input and enables extraction to sequencing in one working day. The clean-up method can be used to generate ultra-long libraries from DNA extracted with phenol chloroform, SDS lysis and CTAB approaches as well as commercial kits such as the NEB Monarch HMW gDNA kits and those from Circulomics (now Pacific Biosciences). Oxford Nanopore Technologies released similar protocols using a spermine-based precipitation buffer (PPT) instead of CoHex and no DNA binding substrate (ONT, ULK001) and then incorporating a star (ONT, ULK114) after our initial protocol release on protocols.io (Cahyani et al. 2021). We compare all these methods here.

Results

The FindingNemo protocol compiles and optimises available ultra-long sequencing protocols

Previously developed protocols

Both the standard and ultra-long (UL) sequencing workflows consist of three main steps: DNA extraction, library preparation, and library clean-up prior to flow cell loading (Figure 1a). Ultra-long sequencing necessitates minimal DNA shearing throughout the library preparation process and during the extraction methods to maintain ultra-high molecular weight (UHMW) DNA. Early experiments pioneering UL sequencing used the transposase-based rapid kit (ONT, e.g., RAD004) to obtain reads with an N50 greater than 100 kb (Jain et al. 2018a). This method eliminated the need for DNA shearing before library preparation. However, the yield from rapid runs was typically low, which might be a consequence of the protocol itself. The lack of purification or clean-up of rapid libraries prior to flow cell loading likely contributed to reduced yields compared to ligation-based sequencing approaches. To better understand the complex interplay of DNA quality, library preparation, read length, and yield, we tested a range of UL protocols aimed at maximising both read length and yield (Figure 1B).

The goal in obtaining UL reads with the rapid kit is to ensure that each DNA molecule is cut only once by a transposase complex, maximising read length. Initially, Quick et al achieved this by saturating the transposase reaction with UHMW DNA (Quick 2018) (Table 1A). An alternative approach taken by Logsdon et al. was to reduce the amount of transposase for a fixed amount of DNA (Logsdon et al. 2021) (Table 1B). However, the sequencing output of this method, when applied to the CHM13 cell line (Logsdon et al. 2021), was lower overall compared to the publicly available nanopore data for GM24385 and GM24631 ligation-based libraries (Supplemental Figure S1A). We hypothesized that this was most likely due to sub-optimal pore occupancy (i.e., the number of available pores sequencing at any time), resulting from the reduced number of adapted DNA ends available for sequencing due to the lower amount of transposase. In testing the Quick protocol, we observed two different read length distributions (see Supplemental Figure S1B, C). Whilst both had N50s of approximately 90 kb, one was dominated by shorter reads. This suggests that the DNA was not cut uniformly. Moreover, although these protocols generated UL reads, they suffered from high variability in output N50, occupancy and yield (Supplemental Figure S1B, C). Typically, as N50 increases, yield and occupancy drop (Supplemental Table S1A, B). The Rocky Mountain protocol (Tyson 2019) (Table 1C) was developed to use the ligation kit (ONT, e.g. LSK109) to address the yield issues experienced with the rapid kit. This protocol utilised varying concentrations of salts and polyethylene glycol (PEG) to precipitate longer reads after the library preparation step. The ligation kit requires light shearing of the input DNA, resulting in non-UL N50s albeit with much improved yields (Tyson 2019). Therefore, this protocol was further modified to use the rapid kit followed by precipitation of the tagmented DNA using a filter paper disc or star as a form of SPRI matrix (KrazyStarFish; Table 1D). The KrazyStarFish (KSF) protocol worked consistently in producing UL N50s, but occupancies and yields were still suboptimal, perhaps a consequence of having DNA clean-up prior to sequencing adapter ligation, consequently producing many free adapters (Supplemental Figure S1D, E; Supplemental Table S1C, D; Methods).

New and developing protocols

ONT released a rapid kit for UL sequencing, SQK-ULK001 (hereafter ULK001) in which the transposase reaction is performed in a large volume (Table 1E), keeping the DNA concentration at

approximately 50 ng/μl. We hypothesised that larger reaction volumes facilitated even diffusion and mixing of DNA and transposase complexes in a more homogeneous state in solution, a feature that we included in our Nemo protocols (Table 1F, G). Efficient transposase reactions result in consistent UL sequencing output as shown by the high rates of occupancy (Supplemental Figure S1F-H, S2, Figure 2, and Supplemental Table S1). Occupancy as a parameter is specifically defined by the percentage of sequencing pores (i.e. pores in 'adapter' and 'strand' states) compared to all pores in available states (i.e. 'adapter', 'strand' and 'pore' state).

We next tested the performance of RAD004 and ULK001 rapid-based kit protocols, modifying parameters such as the number of input cells, the RAD004 dilution buffer, and the clean-up steps whilst maintaining the large volume ratio, and compared these to the control ULK001 protocol (Figure 2A). Regardless of clean-up method (Nanobind disc or glass beads), ULK001-based protocol outputs were comparable in terms of N50s and occupancies (Figure 2B-E). Meanwhile, the output of the library prepared with the RAD004 protocol showed slightly lower performance in general (Figure 2B-E and Supplemental Table S1G).

The total yield of an UL sequencing is not a direct measure of the library quality as yield is also flow cell dependent. The control ULK001 library produced the highest total yield (Figure 2B-C) likely because of the highest pore count of the flow cell (Supplemental Table S1E, Supplemental Figure S3A). Normalising yields between flow cells by taking a ratio of yield to the number of pores utilised during a period of data collection (i.e. yield per pore) allows comparison of run performance. When normalised, the yields per pore of the two ULK001 libraries were comparable (Supplemental Table S1E, F). In all three libraries, the distribution of read lengths against read quality were similar (Supplemental Figure S3B), as well as sequencing speed (Supplemental Figure S3C). These similarities imply that either of the rapid-based kits could be used to successfully produce UL sequencing outputs. Considering all sequencing output parameters, the ULK001 was the optimal protocol for UL sequencing compared to the RAD004 (Supplemental Table S2).

The FindingNemo protocol efficiently purifies and accurately quantifies UHMW DNA and UL libraries

DNA precipitation and recovery: maintaining quality and homogeneity

Prior to the launch of the ULK001 kit we developed a protocol to clean the library after adapter ligation using DNA precipitation chemistry compatible with the motor protein-DNA complex. The filter paper from the KSF protocol was replaced with glass beads, adapting the glass bead's effectiveness to act as a DNA-binding substrate as in the Monarch extraction kit protocol (NEB, T3050). We also surveyed compounds with well-studied DNA precipitating (Ouameur and Tajmir-Riahi 2004; Deng and Bloomfield 1999) and selected spermine and hexamminecobalt(III) chloride (CoHex) for testing DNA precipitation and recovery profiles using glass beads (Supplemental Figure S4A). Optimum final concentrations of the compounds were tested empirically based on previous studies (Pelta et al. 1996).

The ULK001 protocol is carried out in a dilute condition of ~50 ng/μl UHMW DNA. We tested the impact of DNA concentration (10, 20 and 50 ng/μl) on precipitation and recovery rate. CoHex is the best precipitant for both recovery and homogenisation of UHMW DNA at a concentration range of 20 ng/μl to less than 50 ng/μl (Supplemental Figure S4B-D). The optimal DNA recovery of ~75% was obtained using CoHex at DNA concentrations of 20 ng/μl. Spermine shows a significantly higher recovery rate at 50 ng/μl DNA but greater heterogeneity and sample-to-sample variation (Supplemental Figure S4B-D).

The early ULK001 kit version used a spermine-based precipitation buffer (PPT buffer). We tested this alongside our FindingNemo and spermine clean-up protocols, using a DNA input of 40 ng/μl, and split the adapted library into three clean-up preparations. These runs yielded similar UL sequencing metrics (Figure 3A). Albeit with a slightly longer N50, ONT's PPT protocol showed the lowest yield per pore, presumably from the lack of any washing step after the library precipitation that might cause the library to have a higher rate of blocking than the others (Figure 3B).

Homogeneity is measured by calculating the coefficient of variation percentage (%CV) of the DNA concentration. We take at least three concentration measurements from the top, middle and bottom part of the solution using a NanoDrop™ spectrophotometer device. The smaller the %CV, the more

homogeneous, and vice versa. Recovered DNA became more homogeneous after longer in solution, except when precipitated using CoHex at 10 ng/μl DNA concentration (Supplemental Figure S4E). Previous studies showed that a charge neutralisation of 88-90% is required for DNA condensation (Deng and Bloomfield 1999). We hypothesise that at 10 ng/μl DNA, excess CoHex cations affected the water hydration at the DNA-cations interface (Deng and Bloomfield 1999) resulting in a more condensed, difficult-to-dissolve DNA as observed during pipetting. Moreover, NanoDrop is more sensitive when measuring absorbance of the differentially condensed DNA in the solution which explains the larger variation in DNA concentration compared to Qubit™ measurements (Supplemental Figure S4C-D). In contrast, spermine-precipitated DNA shows the least homogeneity in solution (Supplemental Figure S4E). Spermine is also less stable than CoHex; as a solution it is readily oxidised at room temperature (sigmaaldrich.com). CoHex was chosen to precipitate DNA at a concentration range of 20-40 ng/μl UHMW DNA.

The importance of DNA homogeneity and optimum concentration during UL library preparation was shown in the use of a concentrated and non-homogeneous UHMW DNA in sequencing, i.e. more than 60 ng/μl DNA in reaction (Supplemental Figure S5). The sequencing output of this library typifies the output of an undercut library where N50 reaches above 100 kb, but accumulated yield and occupancy are low, as well as a rapid decrease of sequencing pore number (Supplemental Figure S5D). In essence, ensuring the homogeneity of extracted UHMW DNA and setting the reaction concentration below 50 ng/μl provide the optimum conditions for UL sequencing. Longer elution of the library (e.g. overnight) at room temperature with gentle rotation may also help untangle long DNA molecules, facilitating further homogenisation and eventually improving UL sequencing output.

Quantification of UHMW DNA

One of the complications in UL sequencing is in accurately quantitating the UHMW DNA. There is often large variation between the two most used quantification methods: the spectrophotometric and fluorometric based methods, represented by NanoDrop™ and Qubit™ 3 respectively (Supplemental Figure S4F). The viscous nature of UHMW DNA also makes it difficult to withdraw a homogeneous small volume to be measured, even when using a positive displacement pipette. Therefore, we

197 followed a novel Qubit concentration measurement of UHMW DNA to complement the standard
198 NanoDrop method. This method uses Jurkat genomic DNA as the baseline standard for UHMW
199 DNA and a glass bead to homogenise samples before each measurement (Koetsier and Cantor
200 2021). We modified it by combining small volumes from 3-4 different parts of the DNA solution to
201 average out concentration differences within the solution (Cahyani et al. 2021). An accurate
202 concentration measurement requires an ideal ratio of NanoDrop to Qubit DNA values between 1-1.5
203 (Simbolo et al. 2013).

204 As expected, sheared DNA homogenised faster and better than UHMW DNA as shown by the
205 smaller %CV of the sample a few hours after elution (Supplemental Figure S6A). This fits with the
206 fact that shorter DNA molecules will diffuse faster than longer ones. Longer DNA molecules may
207 also have more complex structure and behaviour in solution (Robertson and Smith 2007). The %CV
208 of both types of samples were all below 50% indicating that the FindingNemo method could
209 homogenise DNA solution relatively quickly (Supplemental Figure S6A). It also filtered out
210 contaminants in the input DNA as shown by the increase of the NanoDrop absorbance values after
211 purification (Supplemental Figure S6B) (Koetsier and Cantor 2019). Lastly, the FindingNemo
212 protocol did not interfere with the DNA quality as shown by their size distribution on the pulsed-field
213 electrophoresis gel before and after DNA clean-up (Supplemental Figure S6C).

214 *The FindingNemo protocol vs undercut UL libraries*

215 Combining the FindingNemo approach with either the Quick's or KSF protocols improved
216 sequencing yield but did not significantly improve occupancy or N50 (Supplemental Table S3). We
217 hypothesise that the transposase reaction upstream of the library clean-up is a more critical step in
218 producing optimal ultra-long sequencing outputs. Suboptimal transposase fragmentation may result
219 in long DNA molecules without adapters in the library, due to an insufficient number of transposome
220 complexes compared to DNA molecules, a situation we refer to as an 'undercut'.

221 An undercut library may also occur if the transposome complex-to-DNA ratio is balanced, but the
222 DNA solution is not homogeneous. Viscous, non-homogeneous UHMW DNA is difficult to dilute and
223 can have varying DNA densities at different points in the solution. Consequently, the cutting
224 frequency of the transposome complex will be higher where DNA is more accessible and lower

where long DNA molecules remain in a compacted, less accessible mass. This scenario counterintuitively results in a library with a short N50 despite the high molecular weight nature of the input DNA (Supplemental Figure S8 - 0 rpm sample).

It is also important to consider the interaction between UHMW DNA molecules and the sequencing pores, especially just before a mux scan occurs. During this period, only pores with actively threading (sequencing) DNA molecules will continue sequencing; pores without DNA are inactivated. This process terminates at a fixed time, which is 10 minutes by default. After this period, all pores are inactivated. Molecules still traversing a pore at that point will be ejected, but UHMW DNA increases the risk of stalling or blocking nanopores. At a sequencing speed of 400 bases per second, any molecule longer than 240 kb is at risk of blocking a pore and negatively impacting yield. Therefore, in the FindingNemo protocol, we modified this script parameter ("wind_down") from 600 to 1800 seconds (Cahyani et al. 2024).

The FindingNemo protocol in UHMW DNA extraction

Extractions without kits

We next compared the effect of DNA extraction methods on UL sequencing. The original Quick et al. and Logsdon et al. approaches utilised phenol-chloroform isolation and alcohol-based precipitation with some differences (Table 2A, B). We succeeded in scaling-down these protocols in combination with glass beads to use less cell input and shorten hands-on time to at least half (Table 2C). To reduce the toxicity from phenol use, we also developed phenol-free glass beads-based extraction protocols utilising either SDS, CTAB, or CTAB with CoHex as the lysis buffer (Table 2D, E). These extraction protocols produced (U)HMW DNA, characterised by the distinctly narrow vertical lines/smears and DNA in the wells of the pulsed-field gel electrophoresis results (Supplemental Figure S7).

We also successfully produced UL runs using these DNA samples with varying degrees of occupancies and yields (Figure 4A-J). Only the extraction protocol using a combination of CTAB and CoHex in the lysis buffer could not produce UL N50 (Figure 4F). The inclusion of CoHex during lysis was to test if stable DNA toroid aggregates could be induced (Deng and Bloomfield 1999; Ouameur

and Tajmir-Riahi 2004), so that DNA length could be kept intact during extraction and eventually improve sequencing output. The CoHex precipitated DNA was then washed with ethanol, and it rendered DNA elution from the beads more difficult. Ethanol is required in the wash buffer as it rinses salt, contaminants and other impurities from the cell extract while simultaneously affecting DNA condensation (Arscott et al. 1995; Oda et al. 2016). More optimisation is needed to obtain a consistent rate of elution and recovery in the CTAB-CoHex lysis protocol. Nevertheless, the option of using one of two oppositely charged surfactants as lysis agents, i.e. cationic CTAB (Arseneau et al. 2017) and anionic SDS (Xia et al. 2019), provides flexibility in extracting UHMW DNA from diverse sample types and context.

Extraction using kits

We tested a range of (U)HMW DNA extraction kits during our method development and found that Nanobind CBB Big DNA (Circulomics) and Monarch (NEB) were best suited for our purposes (Table 2F-G). The Nanobind CBB Big DNA Circulomics kit could be used to extract HMW or UHMW DNA, depending on the chosen protocol. DNA extracted with the UHMW protocol showed higher size distribution on a pulsed-field electrophoresis gel compared to the HMW protocol (Supplemental Figure S7A vs S7B). Libraries prepared from UHMW DNA samples extracted using the Circulomics kit produced UL reads and good sequencing metrics (Figure 4I-J).

The Monarch kit protocol initially involved a two-step DNA extraction: first extracting nuclei from cells, then isolating DNA from the nuclei ("Nuclei Prep"). It was later refined to a "Direct Lysis" method, which bypassed the "Nuclei Prep" step by lysing cells and nuclei simultaneously to release DNA. This method yielded high-quality UHMW DNA from GM12878 cells, as shown in pulsed-field gel samples (Supplemental Figure S7A, B). Libraries prepared from these DNA samples indicated that the Direct Lysis approach performed slightly better than the Nuclei Prep method (Supplemental Figure S8).

We also found that the shaking speed during lysis in the Monarch protocol influenced DNA homogeneity, which affected sequencing N50 and occupancy (Supplemental Figure S8; Figure 4G-H; Supplemental Table S4). Lower lysis speeds yielded longer DNA fragments by keeping DNA more intact. However, no shaking during lysis could result in a non-homogeneous UHMW DNA solution, negatively affecting the transposase reaction during library preparation (Supplemental

Figure S8 - 0 rpm). This resulted in a library with non-UL read N50, the lowest occupancy, and yield compared to other libraries (Supplemental Figure S8). The best UL output was achieved by preparing libraries from DNA lysed at a medium speed of 600 to 900 rpm (Supplemental Figure S8; Figure 4G-H; Supplemental Table S4).

Sequencing parameters and performance of the FindingNemo protocol

Maximum yield is obtained at a read N50 between 90-110 kb

On a MinION, yields in excess of 3 Mb per pore after 6 hours of sequencing could be obtained at 93% occupancy, reinforcing that high occupancy is essential for high yield (Figure 4J). However, UL libraries with many short fragments had reduced occupancy as short reads pass quickly through the pores, lowering total yields (Figure 4B, D, E, J). This underscores the importance of homogeneously cut UHMW DNA in UL library preparation. Additionally, targeting an N50 between 90-110 kb is ideal for optimal ultra-long sequencing output, as N50 is inversely proportional to yield (Figure 4J; Supplemental Table S4)

Scalable library loading amount with the FindingNemo protocol

We also tested the impact of the amount of input DNA and library loaded per flow cell. When the input amount of UHMW DNA was less than 5 µg, sequencing output parameters were markedly reduced (Supplemental Table S4, Monarch-01 vs Monarch-05). To produce the optimum N50, occupancy and yield, DNA equivalent to at least one million GM12878 cells (or at least 5 µg) should be used for library preparation, and a minimum of 2 µg of this library loaded onto a MinION flow cell. At least two million cells and 4 µg of library should be used on a PromethION flow cell (Supplemental Table S4, Monarch-01 vs Monarch-02).

FindingNemo in one day: from cells to UL sequencing in one working day

The relatively fast homogenization rate of the FindingNemo clean-up method, combined with the fast DNA extraction when using the Monarch kit can support the preparation of UL libraries from fresh or frozen cells to sequencing in just a working day (Supplemental Figure S9A) (Cahyani 2021). We can obtain UL N50s using this protocol, however yields are markedly reduced (Supplemental Figure S9B-D), likely as a result of the heterogeneity of DNA samples. Yield reductions were exacerbated by the low occupancy level and the abundance of short reads and would be interpreted as undercut libraries (Supplemental Figure S9B-D; Run 2 and 3). Moreover,

the FindingNemo clean-up protocol might not significantly remove these highly abundant short reads due to their high concentration in the libraries.

Using the FindingNemo protocol in diverse sample types, chemistry and platforms

The FindingNemo protocol improved N50, occupancy and yield with ULK114 kit

In 2022, ONT released a new pore and chemistry running on the R10.4.1 flow cells as well as a new ultra-long kit, ULK114. This included new sequencing scripts and a sample rate shift from 4 kHz to 5 kHz. We tested the ULK114 kit according to the manufacturer's instructions and with the FindingNemo protocol on three human cell lines (Figure 5). Libraries prepared with the FindingNemo protocol showed longer N50s, optimal yield, and occupancy (Figure 5A-B). This protocol consistently produced higher yields, regardless of kit chemistry and sequencing programs. Additionally, the 5 kHz script increased the overall mapping rate of reads compared to the 4 kHz program (Figure 5B).

Ratio of transposase to DNA

The ratio of transposase to DNA is crucial for optimising UL sequencing. The ULK001 protocol specifies 1 µl of transposase fragmentation mix (FRA) per 1 million human cells, corresponding to 5-6 µg of UHMW DNA for a 3.2 Gb genome (Table 3). We used this ratio for all GM12878 (Figures 2-5) and HEK293 (Figure 6) runs. Additionally, we tested different FRA to DNA ratios in Mongolian gerbil, tick *Amblyomma variegatum*, and yeast *Saccharomyces cerevisiae* (Table 3). These ratios produced UL sequencing outputs for the gerbil and tick samples (Figure 6A, C). Although the yeast run did not yield UL read N50s (Figure 6B, C), it showed a significant increase in read length compared to existing nanopore-sequenced data (Giordano et al. 2017; Salazar et al. 2017).

Genome ploidy

Genome ploidy anomalies must be considered when preparing UL libraries. For instance, about 4.2% of HEK293 cells are hypotriploid (Synthego HEK293). This means the DNA mass extracted from the same number of cells is higher than from the GM12878 cell line, altering the FRA to DNA ratio (Table 3) and potentially affecting sequencing output by increasing read length and decreasing yield. Sequencing HEK293 DNA with the Monarch kit confirmed this, producing N50s over 190 kb but with lower yields (Figure 6B, C). To maximise yield and N50, it is better to use the

FRA volume to DNA mass ratio rather than the absolute cell number when genome ploidy information is available.

Sample types: to count or to weigh?

The UL protocol requires the same number of human cells per microliter of FRA, whether from a cell line or nucleated blood. However, counting cells in blood samples can be impractical, as seen in the gerbil UL run (Table 3). In such cases, the FRA to DNA mass ratio can be used, provided the DNA is extracted to maintain native chromosome length and properly quantified (Methods). Using this ratio, UL libraries were produced from gerbil blood following the original ULK001 protocol (Figure 6B, C).

The FindingNemo protocol works on different nanopore platforms

We ran pairwise sequencing of identical libraries on both the MinION and PromethION platforms, showing that PromethION produced more yield per pore (Supplemental Table S5). With nearly six times as many pores as a MinION, PromethION is expected to generate more data. The best PromethION results had an N50 of 126 kb and a yield of 3.85 Mb per pore at 88% occupancy (Figure 6C - Tick cells ULK001-Nemo). On the MinION platform, the best results ranged from 106-147 kb N50, with yields of 3.19-3.89 Mb per pore and at least 95% occupancy (Figure 4J). These metrics were all obtained after 6 hours of sequencing.

Relative performance of extraction kits in tick samples

We compared two extraction kits, Circulomics and Monarch, for UL sequencing of a tick cell line with a ~6 Gb genome (Figure 6A, C). Although flow cell occupancies were similar, the Monarch-extracted DNA had longer N50s, higher yields, and a greater proportion of UL reads (Figure 6A right panel, C). The Circulomics library had more short reads, likely due to undercutting (Figure 6C middle panel). The Circulomics-extracted DNA concentration was ~50 ng/μl, while Monarch-extracted DNA was ~26 ng/μl, reflecting differences in sample homogeneity. Moreover, Monarch extractions yielded DNA with a longer size distribution, as shown by pulsed-field electrophoresis gel (Supplemental Figure S10). The Monarch protocol is also more flexible, allowing modifications to target different read-length distributions by adjusting lysis speed (Methods).

Discussion

High-quality, uniformly long DNA is essential for ultra-long reads on Oxford Nanopore sequencers. Our results support the hypothesis that to obtain this length uniformity, transposase reactions should occur in a dilute and homogeneous DNA solution. This can be achieved in a large reaction volume and by thorough mixing to obtain a 'properly cut' library in contrast to 'undercut' as explained before, or 'overcut' which is a higher ratio of FRA to DNA molecules (Figure 7A-B). Additionally, non-adapted DNA molecules may 'crowd' the space around the pores thereby decreasing the diffusion rate of the adapted DNA molecules to be tethered to the pores. The adapted long DNA molecules themselves may inactivate or block the pores during sequencing because of their lengths.

With these challenges, control of DNA quality and quantity is of primary importance. Proper quantification of the number of input cells or tissue weight has direct consequence to controlling DNA yield and to some extent, length distribution. Extended elution and homogenisation are also beneficial; the longer that an UHMW DNA sample is allowed to equilibrate, the more homogeneous it becomes. This reduces the benefit of our FindingNemo in One Day protocol where both maximal read lengths and yields are important. However, as we show, a rapid UL library preparation is possible, albeit with lower total yield. Further optimisation is required, especially at the homogenization of DNA post-extraction, to make this protocol more robust.

Chromosome length distribution of a genome may limit the maximum read N50 that can be obtained, as was the case with the yeast *S. cerevisiae* ultra-long library. Prior knowledge of the genome size and chromosome lengths combined with empirical tests may be required to find the optimum ratio of transposase to DNA amount and produce maximum sequencing N50 and yield. This is a challenging part of ultra-long sequencing as it needs fine tuning and adjustments depending on the sample types used.

We anticipate that the FindingNemo protocol will be generally applicable to tissue samples, could be readily multiplexed, and is likely compatible with automation. Finally, open-source development and sharing of protocols is crucial to enable the widest access to cost-effective high throughput ultra long sequencing and troubleshooting. We established the LongRead Club on protocols.io to house these protocols and promote reproducibility (<https://www.protocols.io/workspaces/long-read->

[club/about](#)). The FindingNemo protocol described here is a valuable addition to the repertoire of ultra-long nanopore sequencing methods available.

Methods

Genomic DNA extraction from cells

Cell sources

GM12878 cells were grown by Darren Crowley (University of Nottingham) and additionally purchased from Coriell Institute for Medical Research, USA. Yeast *S. cerevisiae* cells were given by Stephen Gray (University of Nottingham). Tick *A. variegatum* (AVL/CTVM17) cells were obtained from Alistair Darby's group (University of Liverpool). The HG00233 and HG04054 cells were obtained from Danny Miller's lab (University of Washington, USA). Lastly, HEK293 cells were obtained from New England Biolab (NEB) in Ipswich, MA, USA.

Kit-based extraction

Kits used in the extraction of UHMW DNA were Monarch® HMW DNA Extraction Kit for Cells & Blood (NEB T3050) and Nanobind CBB Big DNA Kit (Circulomics SKU NB-900-001-01) plus the auxiliary Nanobind UL Library Prep Kit (Circulomics SKU NB-900-601-01) according to the manufacturer's protocols with modifications as previously described for the Monarch protocol (Cahyani et al. 2021). The Monarch nuclei prep approach followed the original protocol (NEB T3050) while the direct lysis approach combined the prep and lysis buffers into one step (NEB Direct Lysis). Lysis speed at 600 up to 900 rpm resulted in an optimum UL sequencing output and a minimum of overnight incubation of DNA in solution helped homogenization.

Phenol-chloroform extraction with glass beads

This is a scaled-down version of Quick's original protocol (Quick 2018) and the modifications are as previously described (Cahyani et al. 2021). In summary, pellet of 5 million cells was washed with PBS [Fisher Scientific #15453819] and lysed with an SDS buffer [0.5% SDS, 100mM NaCl, 25mM EDTA pH 8.0, 10mM Tris-HCl pH 8.0] in the presence of RNase A [Qiagen #19101, 20 µg/ml] at 37 °C for 5 min. Proteinase K [Qiagen #19131, 200 µg/ml] was added and incubated at 56 °C for 15 min. Cell lysate was split into two 5PRIME phase-lock gels (Quanta-Bio 2302820). BioUltra TE-saturated phenol (Merck 77607) was added at 1:1 volume ratio. Sample was

homogenised by vertical rotation at 20-30 rpm for 10 min, then centrifuged at 4000 × g for 10 min. The aqueous phase was transferred into new phase-lock gel tubes. DNA purification was repeated with a second round of phenol:chloroform:isoamyl alcohol = 25:24:1 volume ratio. In the third round, only chloroform:isoamyl alcohol=24:1 was used. After the last centrifugation, the aqueous phase from the two tubes was combined in a 5 ml or 15 ml tube and added with 5M ammonium acetate [Sigma-Aldrich #A-7330] (0.4 volume), 3 glass beads (3-mm diameter), and absolute ethanol (2.5 volume). The tube was inverted by hand 20-30 times to precipitate DNA onto the glass beads (or placed in a tube rotator at 10 rpm for 3 min). The supernatant was removed, and bound DNA was washed twice with 60-70% ethanol. Beads with DNA were poured into a bead retainer and quickly spun to remove excess ethanol (or absorbed with a filter paper or tissue). DNA was quickly eluted from beads by pouring them into a 2 ml LoBind [Eppendorf #0030108051] tube containing 10mM Tris-HCl pH 9.0 and incubating at 37 °C for 30 min with regular wide-bore pipette mixing. Incubation was continued overnight at room temperature. Afterwards, using a bead retainer, DNA was spun down at maximum speed for 1 min.

Phenol-free, home brew extractions

This protocol is adapted from the Monarch® HMW DNA Extraction Kit for Cells & Blood as described previously (Cahyani et al. 2021). In summary, a pellet of 1-3 million cells was washed with PBS and lysed in either SDS or CTAB buffer with 200 µg Proteinase K and incubated at 56 °C for 15 min (preferably with shaking at 600-700 rpm). The SDS buffer composition was the same as used in the phenol-based extraction described in the previous section. The CTAB buffer compositions were 2% CTAB, 1.4M NaCl, 25mM EDTA pH 8.0, and 10mM Tris-HCl pH 8.0. After lysis, 100 µg RNase A was added and incubated at 37 °C for 5 min. To precipitate DNA, 5M ammonium acetate was added (0.4 volume), 2-3 glass beads (depending on cell number used), and isopropanol (0.9 volume). Sample was mixed on a rotator at 9 rpm for 5 min (or inverted by hand 20-30 times). Liquid was removed by pipetting and bound DNA was washed twice with 60-70% ethanol. Beads with DNA were poured into a bead retainer and quickly spun to remove excess ethanol (or absorbed with a filter paper or tissue). DNA was quickly eluted with 10mM Tris-HCl pH 8.0 in a 2 ml tube and incubated at 37 °C for 30 min with regular wide-bore pipette mixing.

Incubation was continued overnight at room temperature. Afterwards, using a bead retainer, DNA was spun down at maximum speed for 1 min.

CTAB-CoHex extraction

This protocol is similar to the home-brew CTAB method described above with the following modifications. The CTAB-CoHex buffer compositions used were 2% CTAB, 20mM CoHex, 1.4M NaCl, 25mM EDTA pH 8.0, and 10mM Tris-HCl pH 8.0. To precipitate DNA, sample was diluted with 10mM Tris-HCl pH 8.0, so CoHex concentration became ~5mM and NaCl was ~400mM, *i.e.* ratio of CTAB-CoHex lysis buffer to Tris-HCl was 1:2.5. Then, 2-3 glass beads (depending on cell number used) were added. Sample was mixed on a rotator at 9 rpm for 5 min (or inverted by hand 20-30 times). Liquid was removed by pipetting and bound DNA was washed twice with an ethanol buffer (50% ethanol, 50% of 1mM CoHex). Beads with DNA were poured into a bead retainer and quickly spun to remove excess ethanol (or absorbed with a filter paper or tissue). DNA was quickly eluted with 10mM Tris-HCl pH 8.0 in a 2 ml tube and incubated at 37 °C for 30 min with regular wide-bore pipette mixing. Incubation was continued overnight at room temperature. Afterwards, using a bead retainer, DNA was spun down at maximum speed for 1 min.

Yeast extraction using CTAB and glass beads (Cahyani 2024)

Yeast pellet (~50 µl volume or equivalent to 200 million cells) was washed twice in cold PBS and resuspended in 480 µl SpheroBuffer (1M Sorbitol, 100mM Na₂HPO₄, 100mM EDTA). Lyticase enzyme (5U/µl) was added (20 ul). Sample was incubated at 30°C for 30 minutes and spheroplast formation was checked from time to time. Sample was lysed by adding 480 µl of 2× CTAB lysis buffer (4% CTAB, 2M NaCl, 50 mM EDTA pH 8.0, and 20 mM Tris-HCl pH 9.0), mixed well with wide-bore pipet tip as the buffer is viscous, added with 400 µg Proteinase K (20 mg/ml) and incubated at 56°C for 15 minutes whilst shaking at 700 rpm. To remove RNA, 400 µg of RNase A (100 mg/ml) was added followed by incubation at 56°C for 10 minutes without shaking. Sample was cooled off at room temperature and added with 400 µl of 5M ammonium acetate and 100 µl of 5M NaCl [Thermofisher #AM9759] and mixed by a vertical rotator for 1 minute at 9 rpm. Sample was then centrifuged for 7 minutes at 16,000×g and the lysate supernatant was removed to a new 2 ml tube. Three glass beads and 0.5 ml isopropanol were added, and the sample was mixed on a

vertical rotating mixer at 9 rpm for 5 minutes. Between 0.5-0.7 ml liquid was removed and replaced with the same amount of isopropanol and was rotated for an additional 3 minutes. When DNA had bound tightly around the beads, the liquid was discarded, and bound DNA was washed twice with 1 ml of 70% ethanol. Beads with DNA were poured into a bead retainer and quickly spun to remove excess ethanol (or absorbed with a filter paper or tissue). DNA was quickly eluted with EB (ONT) or 10mM Tris-HCl pH 9.0 in a 2 ml tube and incubated at 37 °C for 30 min with regular wide-bore pipette mixing. Incubation was continued overnight at room temperature. Afterwards, using a bead retainer, DNA was separated from the beads by centrifugation at maximum speed for 1 minute.

Quantification of UHMW DNA

Two nucleic acid quantification methods: fluorometric-based Qubit [Thermo Fisher Scientific, MA, USA] and spectrophotometric-based NanoDrop [Thermo Fisher Scientific], were used in parallel to assess both the quantity and purity of the extracted DNA. The quantification follows the published protocol (Koetsier and Cantor 2021) with some modification as previously described (Cahyani et al. 2021). In brief, DNA is sampled from four different points in the tube: top, upper-middle, lower-middle, and bottom part of the solution. Then, a glass bead was added, and sample was vortexed at full speed for 1 min for concentration measurements. For Qubit measurement, Jurkat genomic DNA (ThermoFisher #SD1111) was used as a concentration standard instead of the lambda DNA provided by the manufacturer (Koetsier and Cantor 2021). To measure DNA homogeneity in solution, %CV of concentration was calculated by directly measuring concentrations of 3-4 different points of the solution using NanoDrop.

Quality check of UHMW DNA using PFGE

To assess DNA quality and read-length distribution, 200-300 ng extracted DNA was run on 0.75% gold agarose gel (SeaKem #50150). Pulsed-field gel electrophoresis was carried out using Pippin Pulse device (Sage Science) at a pre-set mode of “5-430 kb” for 16 hours. Gel was stained with RedSafe DNA stain (ChemBio #21141) and visualised with a Bio-Rad GelDoc imaging system.

Rapid kit-based UL library preparation

Ultra-long kit-Nemo

The ONT ultra-long kits, i.e. ULK001 and ULK114, are based on the cleaving and tagging of the DNA by the transposome complexes to insert transposase adapters, which then are ligated to the sequencing adapters (ONT ULK001). Our modifications of the ULK001 method were described in previously published protocol (Cahyani et al. 2021). To summarise, DNA extracted from 1-6 million cells was diluted to a final concentration of 20-40 ng/μl (in the final reaction volume). All mixing was done with wide-bore pipette tips. Meanwhile on ice, the dilution buffer (FDB) was well mixed with FRA (1 μl per 1 million human cells DNA), then combined to the DNA solution, and quickly incubated at 23 °C for 5-10 min (depending on the reaction volume). Enzyme was inactivated at 70 °C for 5 min and cooled down to room temperature. Sequencing adapter (RAP-F) was added (0.83 μl per 1 μl FRA used) and the sample was incubated at room temperature for at least 30 min. The ULK114 modifications are similar to ULK001 and are described in the newly published and updated protocol also on protocols.io (Cahyani et al. 2024). In summary, library preparation to run on a PromethION flow cell needs extracted DNA from 6 million cells that was diluted to a final concentration of 20-40 ng/μl (in the final reaction volume). On ice, the dilution buffer (FDB) was well mixed with FRA (6 μl per 6 million human cells DNA), then combined to the DNA solution, and quickly incubated at 23 °C for 10 min with 9 rpm rotation. Enzyme was inactivated at 75 °C for 10 min and cooled down to room temperature. Sequencing adapter (RA) was added (5 μl per 6 μl FRA) and the sample was rotated at 9 rpm for 30 minutes.

RAD004-Nemo

The RAD004-Nemo protocol followed the same route as the ULK001-Nemo protocol described previously (Cahyani et al. 2021). However, the transposase and adapter used were the ones from the RAD004 kit (*i.e.* FRA and RAP, respectively). The dilution buffer used to replace FDB was the 4× MuA buffer [100 mM Tris-HCl pH 8.0, 40 mM MgCl₂, 440 mM NaCl, 0.2% Triton X-100, 40% Glycerol].

Krazy StarFish (KSF)

The KSF protocol is an alternative UL sequencing that uses a filter paper instead of glass beads. It is a cheap, quick-and-dirty method to produce UL sequencing based on the RAD004 kit. This

method was described previously (Cahyani et al. 2021) and summarised as follows. An amount of ~7.5 µg DNA in 75 µl volume was mixed well with 25 µl 4× MuA buffer. Meanwhile in another tube on ice, 25 µl 4× MuA buffer was diluted with 74 µl water and added with 1 µl FRA (or up to 1.2 µl FRA if DNA was quite viscous). These two solutions were mixed well on ice, aliquoted into two PCR tubes, and treated at 30° C for 1 min, 80° C for 1 min then cooled to room temperature. Reactions were pooled into a single tube, 12 µl 5M NaCl added, and mixed gently. A Krazy StarFish filter (made from Whatman filter paper no. 3 with a star hole puncher) was submerged in the solution, 142 µl isopropanol added, and mixed by inversion 20-30 times allowing UHMW DNA to collect and condense onto the filter. Liquid was removed and filter paper was washed twice with 60-70% ethanol. Excess ethanol was removed by pipetting and paper was air-dried for 2 min. After the paper was moved to a clean tube, the library was eluted by adding 120 µl of elution buffer [10mM Tris-HCl pH 8.0]. Before sequencing, 37.5 µl library was added with 0.5 µl RAP and incubated for at least 30 min at room temperature.

Nemo Clean Up: UHMW DNA Library Purification with Glass Beads and CoHex

The Nemo protocol is at the heart of this study and was previously described as an alcohol-free library clean-up method (Cahyani et al. 2021, 2024). Typically, for DNA equivalent to 1 million cells or more, 3 glass beads were used. In brief, glass beads were added to the library to be cleaned-up in a 1:1 volume ratio of 10 mM CoHex (2×), i.e. final CoHex concentration is 5 mM. Sample was rotated at 9 rpm for 3 min, or 20-25 times if using hand inversion. Liquid was removed and beads were washed twice with PEGW buffer [10% PEG-8000, 0.5M NaCl] for 3 minutes each. Excess buffer was removed with a fine pipette tip after a quick spin of the tube. Alternatively, the beads with DNA could be poured into a bead retainer (home-made or using the NEB-T3004) and the buffer was absorbed by tapping on a filter paper/tissue. Beads should not be let dry, were quickly eluted with 10mM Tris-HCl pH 9.0 or EB (ONT) and incubated at 37 °C for 30 min with regular wide-bore pipette mixing. Incubation was continued overnight (preferable) at room temperature. Afterwards, using a bead retainer on a new tube, DNA was spun down at maximum speed for 1 min.

Flowcell Priming and Library Loading

All flow cells were primed according to the ULK001 and ULK114 protocols. Library was loaded onto a primed FLO-MIN106 or FLO-PRO002 (R9.4.1), or FLO-MIN114 or FLO-PRO114 (R10.4.1) flow cells for sequencing on a GridION or PromethION, respectively, using the accompanying MinKNOW version software. Loaded UL library was let to tether on the flow cell for at least 30 min before run was started using the ONT ultra-long script with mux scan of every 6 hours for R9.4.1 and the default 1.5 hours for R10.4.1.

FindingNemo in one day

To run UL sequencing within the same day as the DNA extraction, FindingNemo *in one day* protocol was developed as previously described (Cahyani 2021). In brief, cells were extracted using the Monarch kit and lysed at 700-800 rpm. Eluted DNA was incubated at 37 °C and regularly mixed for at least two hours. After quantification showed that DNA was sufficiently homogeneous (%CV < 100%), the library was prepared following the ULK001-Nemo protocol. Eluted library was further incubated at 37 °C for at least 30 min with regular pipette mixing, quantified and loaded on a primed flow cell.

NanoPlot Data Analysis and Plotting

Violin plots of sequencing speed and distribution graphs comparing read-lengths and read quality were created using the summary file data processed by the NanoPlot tool (De Coster et al. 2018).

R Data Analyses and Plotting

Read-length histograms, occupancy, yield and other sequencing metrics were analysed from the summary, mux, and duty time (or in newer filename version, pore activity) output files using in-house R (R Core Team 2024) scripts available on GitHub (<https://github.com/cinswasti/FindingNemo>) or in the Supplemental Codes.

Sequence sources and alignment

For human work, the reference genome used is the hg38 (NCBI RefSeq assembly GCF_000001405.40). Sequence data were aligned to this reference genome using minimap2 (Li 2021). Nanopore data of GM24385 and GM24631 are publicly available on AWS registry and can

be accessed using the command line `aws s3 ls --no-sign-request s3://ont-open-data`.`

Data Access

The FASTQ data generated in this study had been submitted to the European Nucleotide Archive (ENA; <https://www.ebi.ac.uk/ena/browser/home>) under accession number PRJEB76809. The FindingNemo protocol can be found in Supplemental_Protocol_FindingNemo or on the protocols.io website (<https://dx.doi.org/10.17504/protocols.io.5jyl8p38rg2w/v1>).

Competing Interest Statement

M.L., J.T, J.Q. and N.L were members of the MinION early access program and have received free flow cells and sequencing reagents in the past. All have received reimbursement for travel and accommodation to speak at events organised by ONT.

Acknowledgements

We would like to thank Giron Koetsier (NEB), Simon Mayes (ONT), and Kelvin Liu (Circulomics) for lending their expertise and/or advanced product trials. We thank Darren Crowley and Luke Simpson for supplying the GM12878 cells, Stephen Gray for the yeast cells, Danny Miller and Miranda Galey (University of Washington) for the Human Genome (HG) cells, Lesley Bell-Sakyi and Alistair Darby (University of Liverpool) for the tick *A. variegatum* cells, Giron Koetsier (NEB) for the HEK293 cells, and John Mulley (Bangor University) for the use of the gerbil's sequencing metrics. John Tyson is grateful to Terrance Snutch. This work was supported by the Wellcome Trust (grant number 204843/Z/16/Z, I.C., N.H., N.L and M.L.).

Author contributions

N.L, M.L, J.T and J.Q conceived the project. I.C. led the work with contributions from J.Q, J.T, N.H and C.M. I.C drafted the manuscript with input from all authors. All authors reviewed, revised and approved the final version of the manuscript.

References

- Allers T, Lichten M. 2000. A method for preparing genomic DNA that restrains branch migration of Holliday junctions. *Nucleic Acids Res* **28**: 6.
- Arscott PG, Ma C, Wenner JR, Bloomfield VA. 1995. DNA condensation by cobalt hexaammine(III) in alcohol–water mixtures: Dielectric constant and other solvent effects. *Biopolymers* **36**: 345–364. <https://onlinelibrary.wiley.com/doi/10.1002/bip.360360309> (Accessed February 25, 2022).
- Arseneau J-R, Steeves R, Laflamme M. 2017. Modified low-salt CTAB extraction of high-quality DNA from contaminant-rich tissues. *Mol Ecol Resour* **17**: 686–693. <https://onlinelibrary.wiley.com/doi/10.1111/1755-0998.12616>.
- Cahyani I. 2021. FindingNemo in OneDay: Ultra-Long ONT Library Preparation from Cell to Flowcell in One Day. <https://www.protocols.io/view/findingnemo-in-oneday-ultra-long-ont-library-prepa-14egnzqzyg5d/v1> (Accessed November 4, 2024).
- Cahyani I, Tyson J, Holmes N, Quick J, Loman N, Loose M. 2021. FindingNemo: A Toolkit of CoHex- and Glass Bead-based Protocols for Ultra-Long Sequencing on ONT Platforms. *protocols.io*. <https://www.protocols.io/view/findingnemo-a-toolkit-of-cohex-and-glass-bead-base-dm6gpwkn1lzp/v1> (Accessed March 11, 2024).
- Cahyani I, Tyson J, Holmes N, Quick J, Loman N, Loose M. 2024. FindingNemo (v.kit14): A Toolkit for DNA Extraction, Library Preparation and Purification for Ultra Long Nanopore Sequencing. *protocols.io*. <https://www.protocols.io/view/findingnemo-v-kit14-a-toolkit-for-dna-extraction-l-5jyl8p38rg2w/v1/metadata> (Accessed April 17, 2024).
- De Coster W, D’Hert S, Schultz DT, Cruts M, Van Broeckhoven C. 2018. NanoPack: visualizing and processing long-read sequencing data. *Bioinformatics* **34**: 2666–2669. <https://dx.doi.org/10.1093/bioinformatics/bty149> (Accessed November 4, 2024).
- Deng H, Bloomfield VA. 1999. Structural Effects of Cobalt-Amine Compounds on DNA Condensation. *Biophys J* **77**: 1556–1561. <https://linkinghub.elsevier.com/retrieve/pii/S0006349599770037> (Accessed March 11, 2024).

649 Giordano F, Aigrain L, Quail MA, Coupland P, Bonfield JK, Davies RM, Tischler G, Jackson
650 DK, Keane TM, Li J, et al. 2017. De novo yeast genome assemblies from MinION,
651 PacBio and MiSeq platforms. *Sci Rep* **7**: 1–10.

652 Jain M, Koren S, Miga KH, Quick J, Rand AC, Sasani TA, Tyson JR, Beggs AD, Diltthey AT,
653 Fiddes IT, et al. 2018a. Nanopore sequencing and assembly of a human genome with
654 ultra-long reads. *Nat Biotechnol* **36**: 338–345. <http://www.nature.com/articles/nbt.4060>.

655 Jain M, Olsen HE, Turner DJ, Stoddart D, Bulazel K V, Paten B, Haussler D, Willard HF,
656 Akesson M, Miga KH. 2018b. Linear assembly of a human centromere on the Y
657 chromosome. *Nat Biotechnol* **36**: 321–323. <https://www.nature.com/articles/nbt.4109>
658 (Accessed March 11, 2024).

659 Kankia BI, Buckin V, Bloomfield VA. 2001. Hexaminecobalt(III)-induced condensation of
660 calf thymus DNA: circular dichroism and hydration measurements. *Nucleic Acids Res* **29**:
661 2795. [/pmc/articles/PMC55774/](https://pubmed.ncbi.nlm.nih.gov/11611111/) (Accessed July 19, 2022).

662 Koetsier G, Cantor E. 2019. A Practical Guide to Analyzing Nucleic Acid Concentration and
663 Purity with Microvolume Spectrophotometers. *New England Biolabs* **1**: 1–8.
664 [https://www.neb.com/-/media/catalog/application-](https://www.neb.com/-/media/catalog/application-notes/mvs_analysis_of_na_concentration_and_purity.pdf?rev=be7c8e19f4d34e558527496ea51623dc)
665 [notes/mvs_analysis_of_na_concentration_and_purity.pdf?rev=be7c8e19f4d34e5585274](https://www.neb.com/-/media/catalog/application-notes/mvs_analysis_of_na_concentration_and_purity.pdf?rev=be7c8e19f4d34e558527496ea51623dc)
666 [96ea51623dc](https://www.neb.com/-/media/catalog/application-notes/mvs_analysis_of_na_concentration_and_purity.pdf?rev=be7c8e19f4d34e558527496ea51623dc).

667 Koetsier PA ‘Giron’, Cantor EJ. 2021. A simple approach for effective shearing and reliable
668 concentration measurement of ultra-high-molecular-weight DNA. *Biotechniques* **71**: 439–
669 444. <https://www.future-science.com/doi/10.2144/btn-2021-0051> (Accessed August 4,
670 2021).

671 Logsdon G. 2020. HMW gDNA purification and ONT ultra-long-read data generation v1.
672 [https://www.protocols.io/view/hmw-gdna-purification-and-ont-ultra-long-read-data-](https://www.protocols.io/view/hmw-gdna-purification-and-ont-ultra-long-read-data-bchhit36)
673 [bchhit36](https://www.protocols.io/view/hmw-gdna-purification-and-ont-ultra-long-read-data-bchhit36) (Accessed December 13, 2024).

674 Logsdon GA, Vollger MR, Hsieh PH, Mao Y, Liskovych MA, Koren S, Nurk S, Mercuri L,
675 Dishuck PC, Rhie A, et al. 2021. The structure, function and evolution of a complete
676 human chromosome 8. *Nature*. <http://dx.doi.org/10.1038/s41586-021-03420-7>.

- Oda Y, Sadakane K, Yoshikawa Y, Imanaka T, Takiguchi K, Hayashi M, Kenmotsu T, Yoshikawa K. 2016. Highly Concentrated Ethanol Solutions: Good Solvents for DNA as Revealed by Single-Molecule Observation. *Chemphyschem* **17**: 471. [/pmc/articles/PMC4770436/](https://pubmed.ncbi.nlm.nih.gov/2670436/) (Accessed February 25, 2022).
- Ouameur AA, Tajmir-Riahi H-A. 2004. Structural Analysis of DNA Interactions with Biogenic Polyamines and Cobalt(III)hexamine Studied by Fourier Transform Infrared and Capillary Electrophoresis. *Journal of Biological Chemistry* **279**: 42041–42054. <https://linkinghub.elsevier.com/retrieve/pii/S002192582072675X> (Accessed February 25, 2022).
- Pelta J, Livolant F, Sikorav J-L. 1996. DNA Aggregation Induced by Polyamines and Cobalthexamine. *Journal of Biological Chemistry* **271**: 5656–5662. <https://linkinghub.elsevier.com/retrieve/pii/S0021925818979068> (Accessed March 11, 2024).
- Quick J. 2018. Ultra-long read sequencing protocol for RAD004. *Protocols* **10**: 1–16. <https://www.protocols.io/view/ultra-long-read-sequencing-protocol-for-rad004-mrvc566> (Accessed August 4, 2021).
- R Core Team (2024). R: A language and environment for statistical computing. R Foundation for Statistical Computing, Vienna, Austria. URL: <https://www.R-project.org/>.
- Robertson RM, Smith DE. 2007. Strong effects of molecular topology on diffusion of entangled DNA molecules. *Proc Natl Acad Sci U S A* **104**: 4824–4827. www.pnas.org/cgidoi/10.1073/pnas.0700137104 (Accessed March 12, 2024).
- Rodrigue S, Materna AC, Timberlake SC, Blackburn MC, Malmstrom RR, Alm EJ, Chisholm SW. 2010. Unlocking Short Read Sequencing for Metagenomics. *PLoS One* **5**: e11840. <https://journals.plos.org/plosone/article?id=10.1371/journal.pone.0011840> (Accessed February 7, 2022).
- Salazar AN, Gorter de Vries AR, van den Broek M, Wijsman M, de la Torre Cortés P, Brickwedde A, Brouwers N, Daran J-MG, Abeel T. 2017. Nanopore sequencing enables near-complete de novo assembly of *Saccharomyces cerevisiae* reference strain

CEN.PK113-7D. *FEMS Yeast Res* **17**.

<https://academic.oup.com/femsyr/article/17/7/fox074/4157789> (Accessed July 22, 2021).

Simbolo M, Gottardi M, Corbo V, Fassan M, Mafficini A. 2013. DNA Qualification Workflow for Next Generation Sequencing of Histopathological Samples. *PLoS One* **8**: 62692.

www.plosone.org (Accessed March 12, 2024).

Stortchevoi A, Kamelamela N, Levine SS. 2020. SPRI beads-based size selection in the range of 2-10kb. *Journal of Biomolecular Techniques* **31**: 7–10.

[/pmc/articles/PMC6944320/](https://pubmed.ncbi.nlm.nih.gov/3244320/) (Accessed February 7, 2022).

Synthego HEK293. HEK293 Cells: Background, Applications, Protocols, and More.

[synthego.com](https://www.synthego.com). <https://www.synthego.com/hek293> (Accessed March 12, 2024).

Tyson J. 2019. Rocky Mountain adventures in Genomic DNA sample preparation, ligation protocol optimisation / simplification and Ultra long read generation v1.

<https://www.protocols.io/view/rocky-mountain-adventures-in-genomic-dna-sample-pr-7euhjew> (Accessed November 4, 2024).

Xia Y, Chen F, Du Y, Liu C, Bu G, Xin Y, Liu B. 2019. A modified SDS-based DNA extraction method from raw soybean. *Biosci Rep* **39**: 1–10.

<https://portlandpress.com/bioscirep/article/39/2/BSR20182271/110910/A-modified-SDS-based-DNA-extraction-method-from>.

Figure legends

Figure 1. Summary of the FindingNemo toolkit: (A) Similarities and differences between the standard and ultra-long (UL) protocol that centre on the way UHMW DNA is handled. Created with BioRender (<https://biorender.com>). (B) The route of FindingNemo toolkit consists of five main steps, with options at the extraction, library preparation and clean-up steps. Options with asterisks (*) are our current lab workflow for the ultra-long protocol.

Figure 2. Validation of the FindingNemo protocol - sequencing outputs between rapid-based library preparation protocols (ULK001-Nemo and RAD004-Nemo protocol) showing comparable

performance when compared to the control ULK001 protocol. (A) Workflow of the library preparation protocols, (B) Total yields sub categorised by read lengths, (C) Time lapse of yields, (D) Read length distributions; dashed black vertical lines denote N50s, (E) Time lapse of relative pore occupancies. Each library was loaded three times (i.e. total DNA equivalent to 3 million GM12878 cells) on a MinION R.9.4.1 flow cell following a nuclease flush protocol. All data are after 65 hours of sequencing on the GridION platform.

Figure 3. Choosing the best precipitating agent that ensures fast homogenisation: (A) Read length distribution of libraries precipitated with different buffers; dashed black vertical lines denote N50s, (B) Sequencing metrics after 6 hours of run (excluding the first 10 minutes). All libraries were extracted from GM12878 cells using Monarch direct lysis protocol with an input DNA concentration of ~40 ng/μl split into three clean-up preparations, loaded on MinION R.9.4.1 flow cells and sequenced on GridION platform.

Figure 4. Combinatorial effects of extraction and library clean-up methods, either with or without the use of kits, on the UL sequencing outputs: (A-I) Read length distributions of the libraries: graph titles denote the used extraction kit/method and the clean-up protocol, dashed black vertical lines denote N50s. (J) Sequencing metrics after 6 hours of run (excluding the first 10 minutes). DNA was extracted from GM12878 cells. Each library was loaded on a MinION R9.4.1 flow cell and run on the GridION platform.


Figure 5. Ultra-long N50s were obtained using the new ULK114 kit from three cell lines: (A) Read length distribution of GM12878 library run with the earlier ULK114 script version (left panel), while the other two libraries (middle and right panel) were run with the most current 5 kHz version of the script. Vertical dashed black lines indicate N50s. (B) Run metrics of the three libraries showing the Nemo protocol performed better than ONT's.

Figure 6. Ultra-long sequencing output and metrics of non-GM12878 cells and different sequencing platforms: gerbil blood, tick *Amblyomma variegatum*, yeast *Saccharomyces cerevisiae*


and human HEK293 cells using either the ULK001 or the ULK001-Nemo protocol (as labelled). (A) Read length distributions of libraries run on PromethION R9.4.1 flow cells; dashed black vertical lines denote N50s, (B) Read length distribution of the HEK293 and yeast *S. cerevisiae* libraries run on MinION R9.4.1 flow cells, (C) Sequencing metrics of the libraries. Each metric was shown after 6 hours of run (excluding the first 10 minutes).






Figure 7. Optimum conditions for ultra-long sequencing - to obtain an optimum UL sequencing output, there needs to be a 1:1 ratio of FRA to DNA molecules in a homogeneous solution ascertained by proper cell counting and/or DNA mass measurement. Image was created with BioRender (<https://biorender.com>).

A

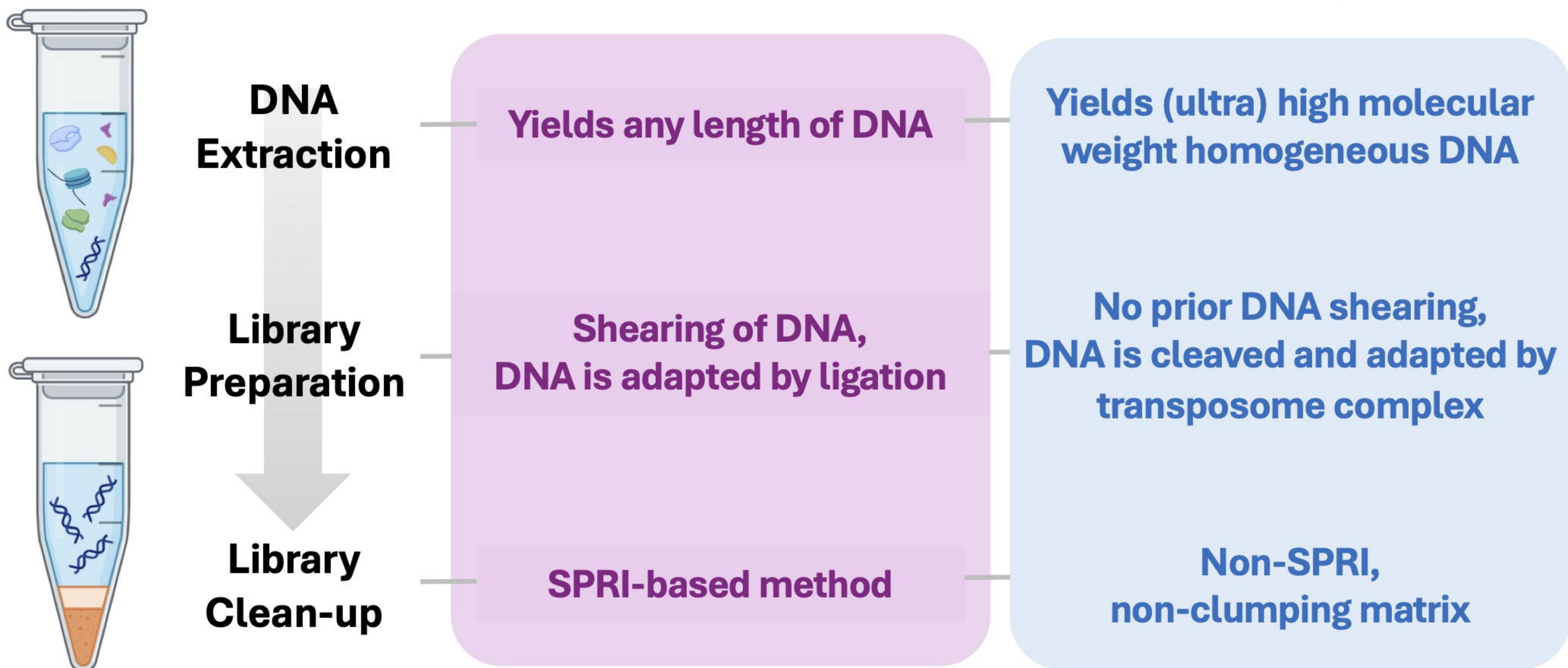


B

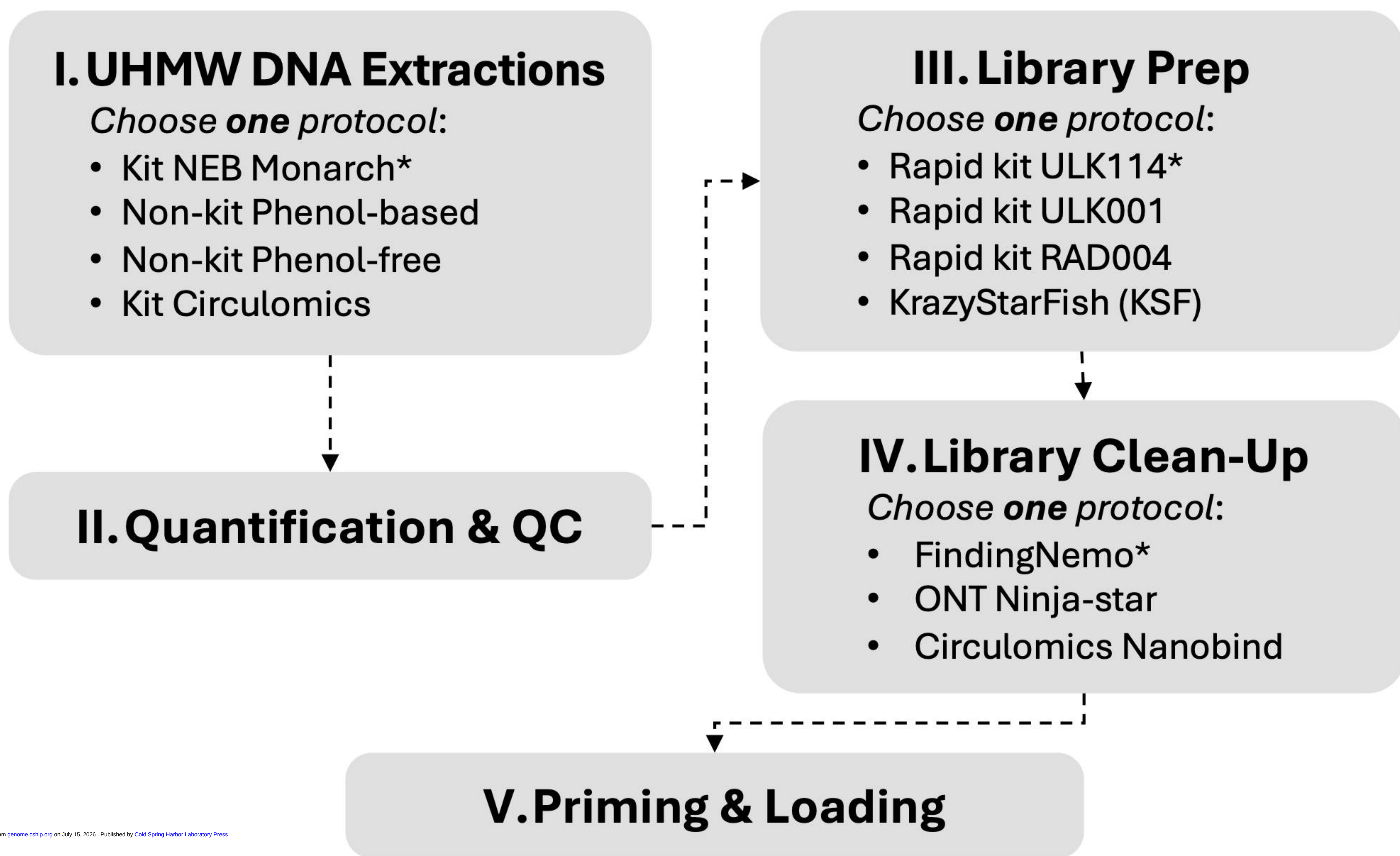

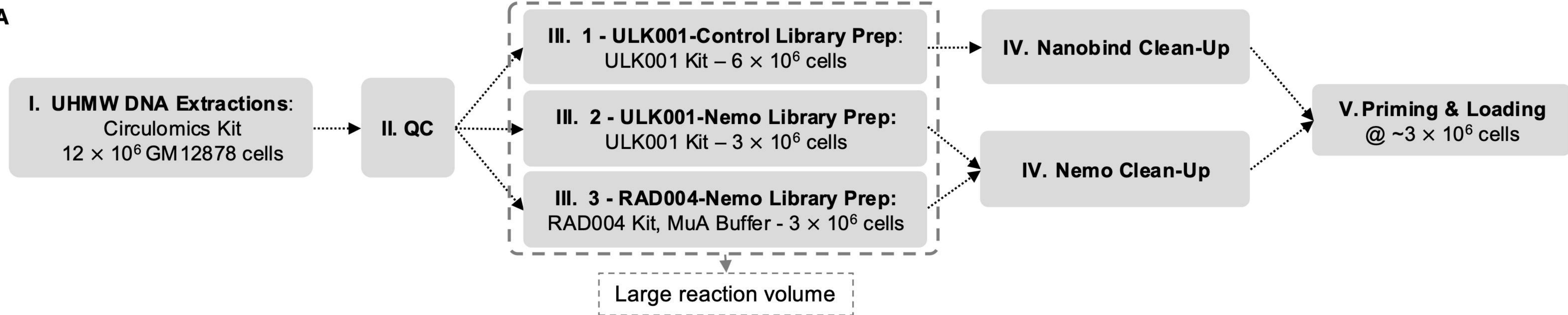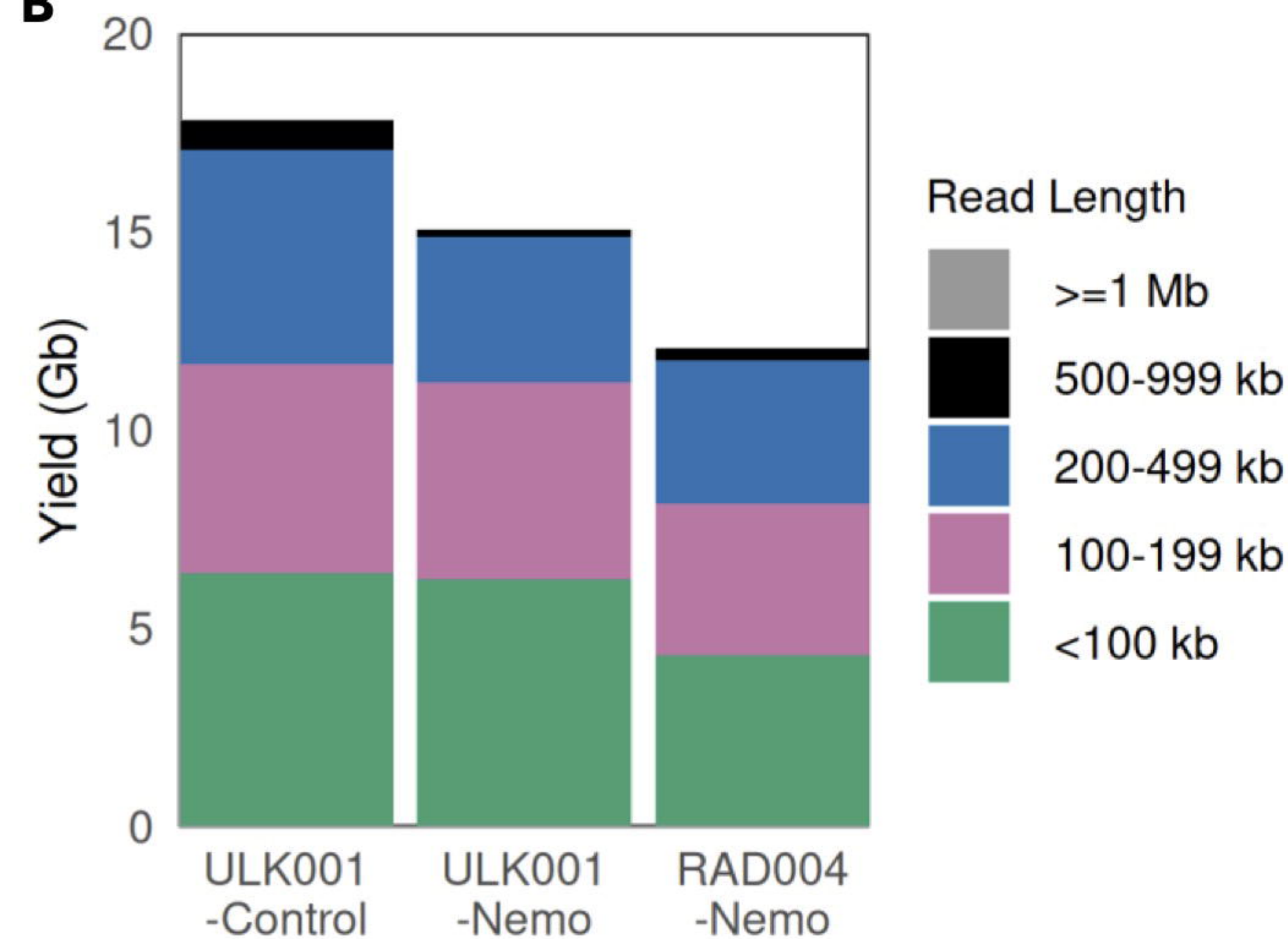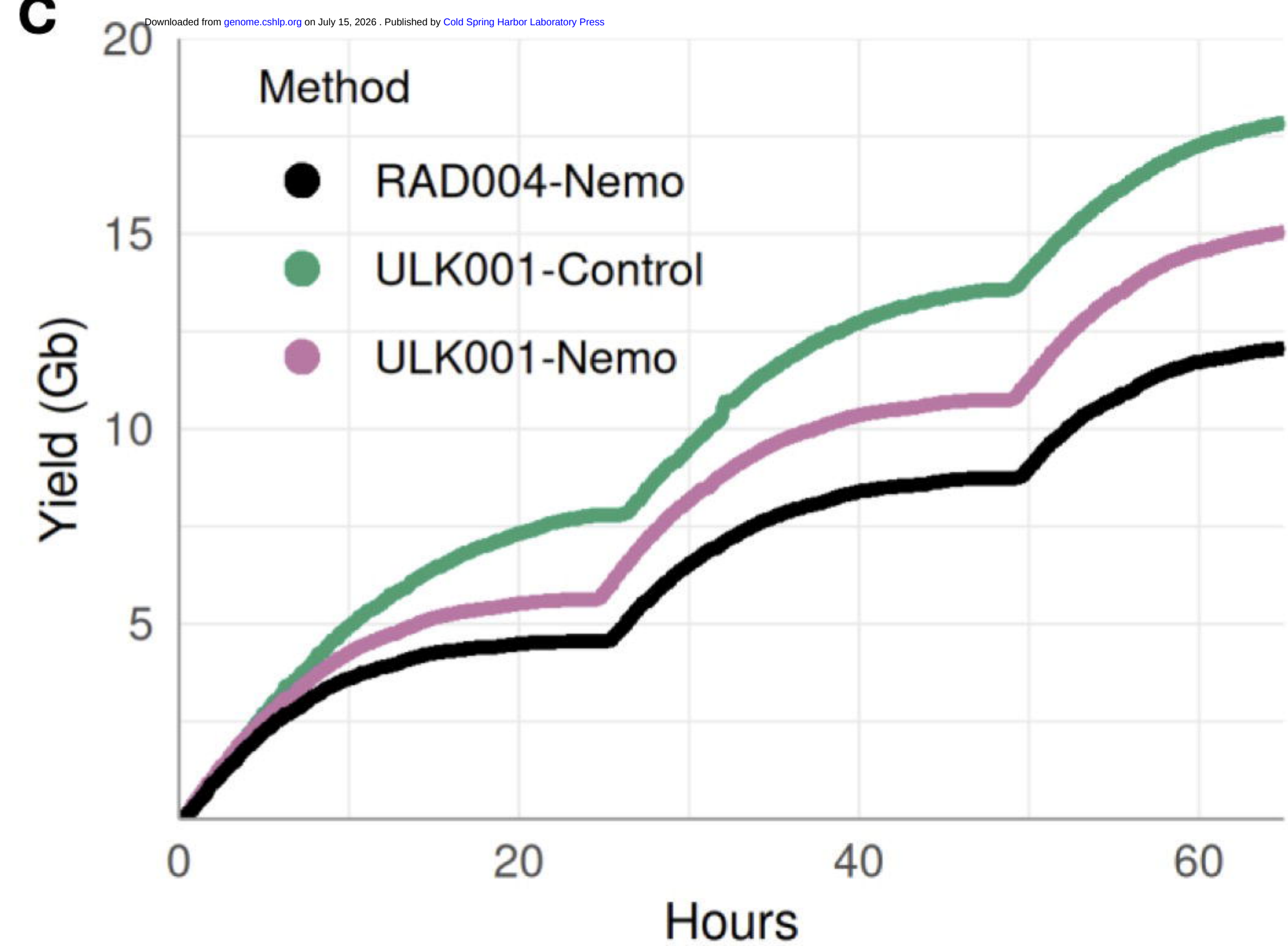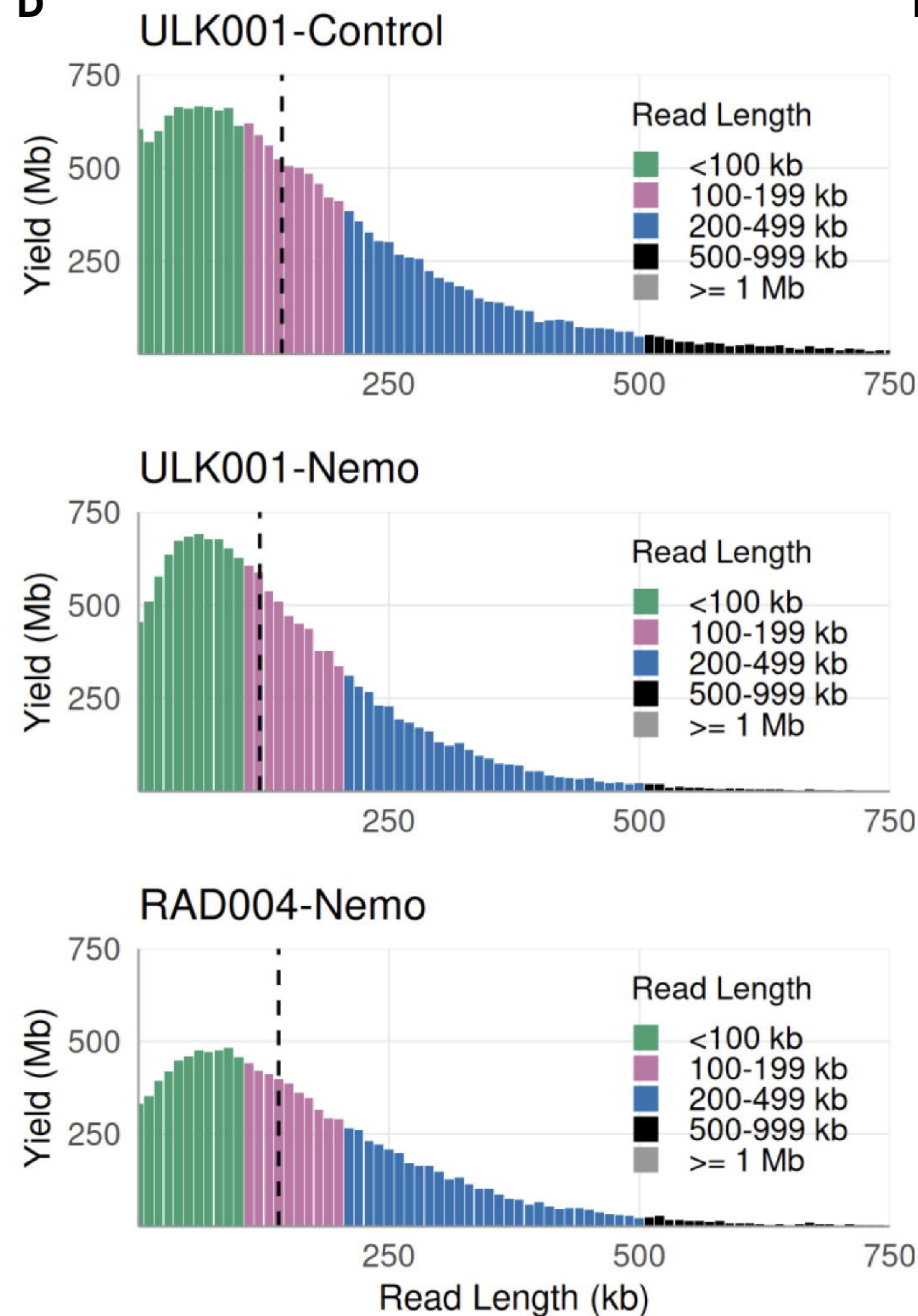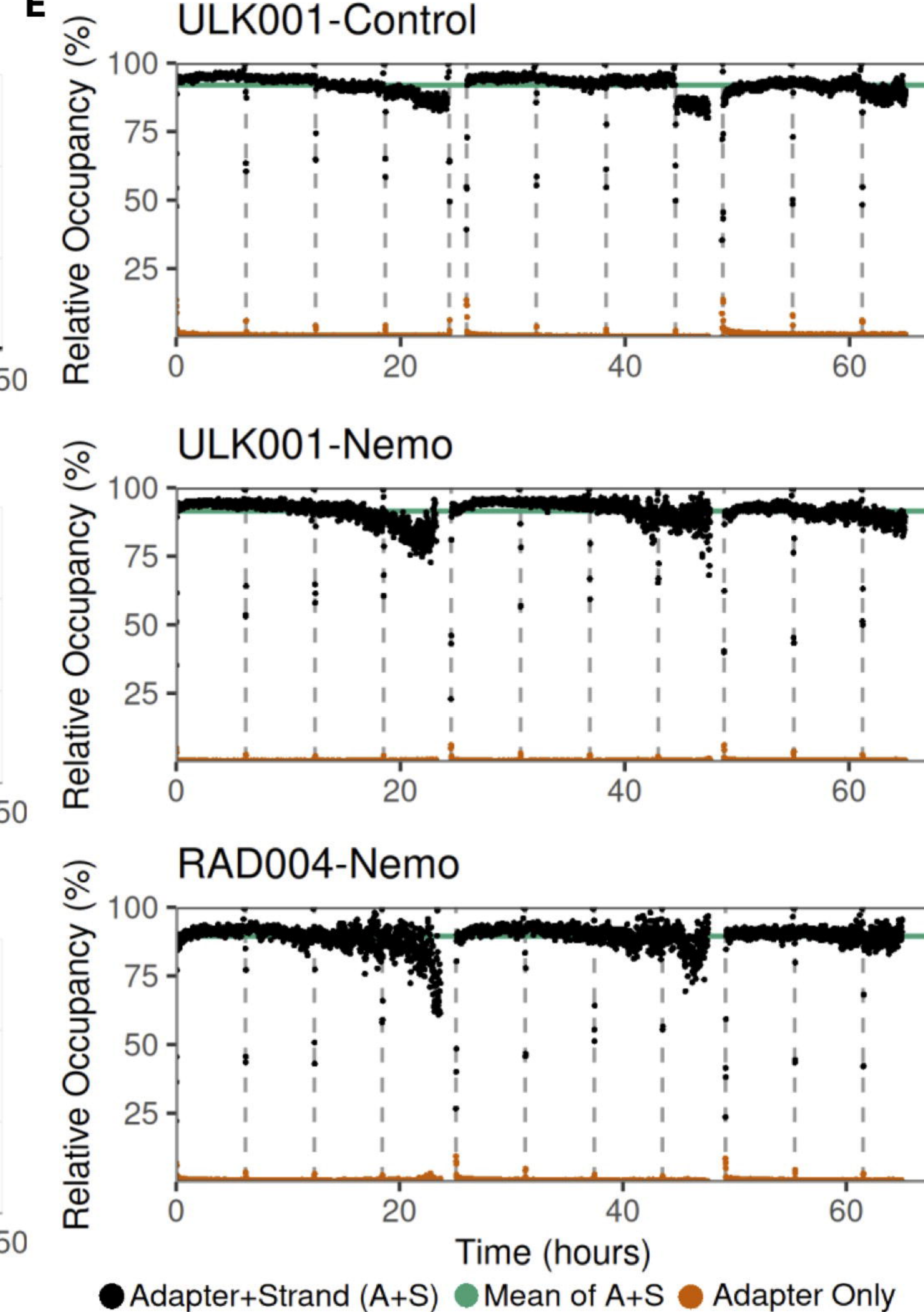A**B****C****D****E**

A

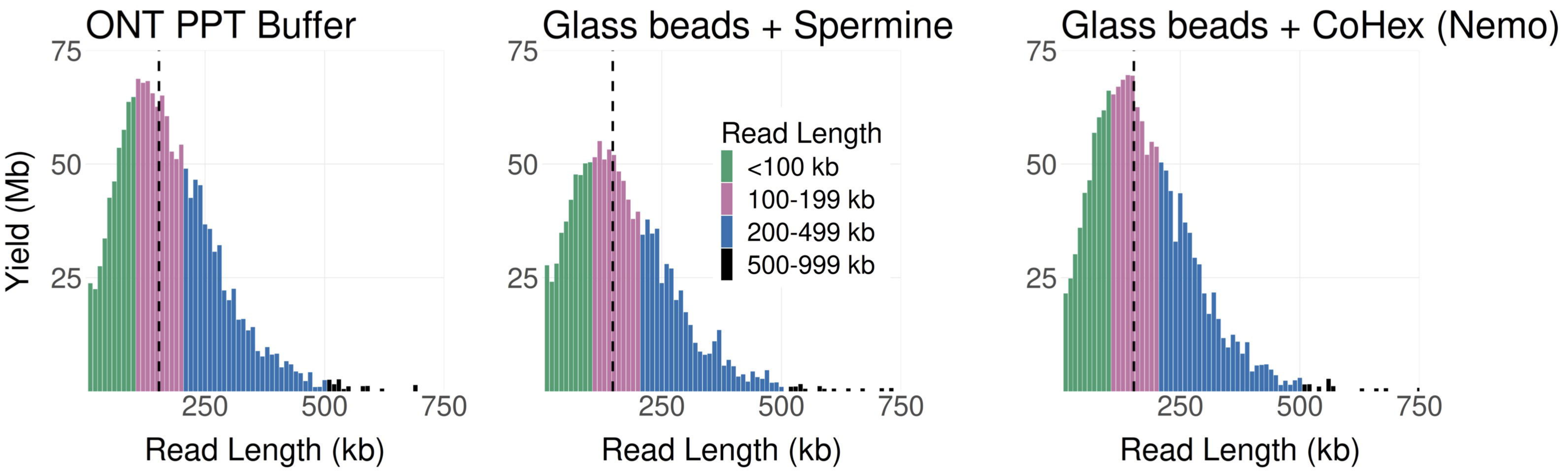


B


Sequencing Metrics at 6 Hours of Run

Downloaded from genome.cshlp.org on July 15, 2025. Published by Cold Spring Harbor Laboratory Press

Sample	N50 (bp)	Pores Used	Total Yield (% Bases $\geq$ 100 kb)	Yield per Pore (Mb)	Average Occupancy (%)	% CV of Library Concentration (Nanodrop)
ONT PPT buffer	153,295	1010	1.59 Gb (72.55%)	1.57	85.65	14.8
Glass Beads + Spermine	147,386	635	1.28 Gb (69.51%)	2.01	85.32	11.1
Glass Beads + CoHex (Nemo)	152,653	794	1.61 Gb (72.20%)	2.03	89.29	1.3

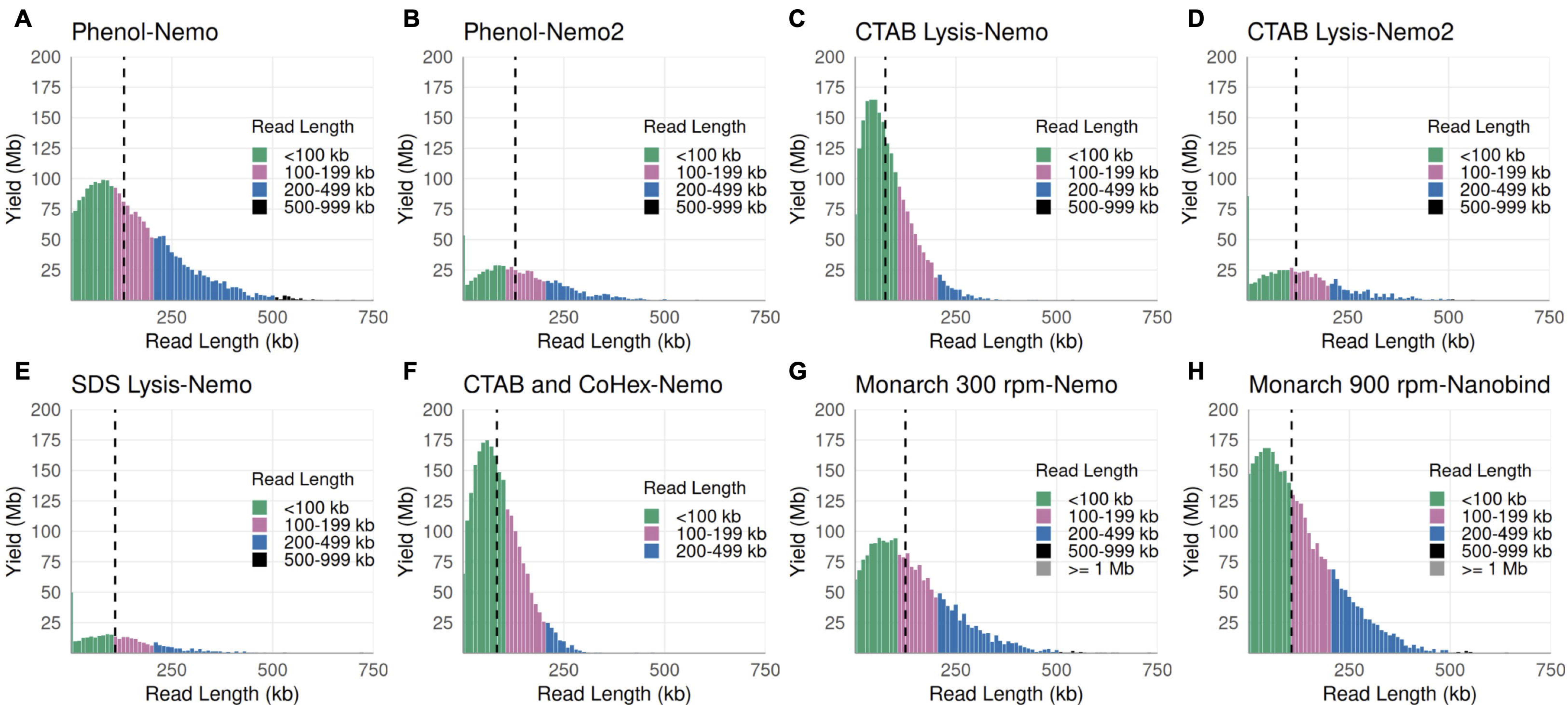


Downloaded from genome.cshlp.org on July 15, 2025. Published by Cold Spring Harbor Laboratory Press

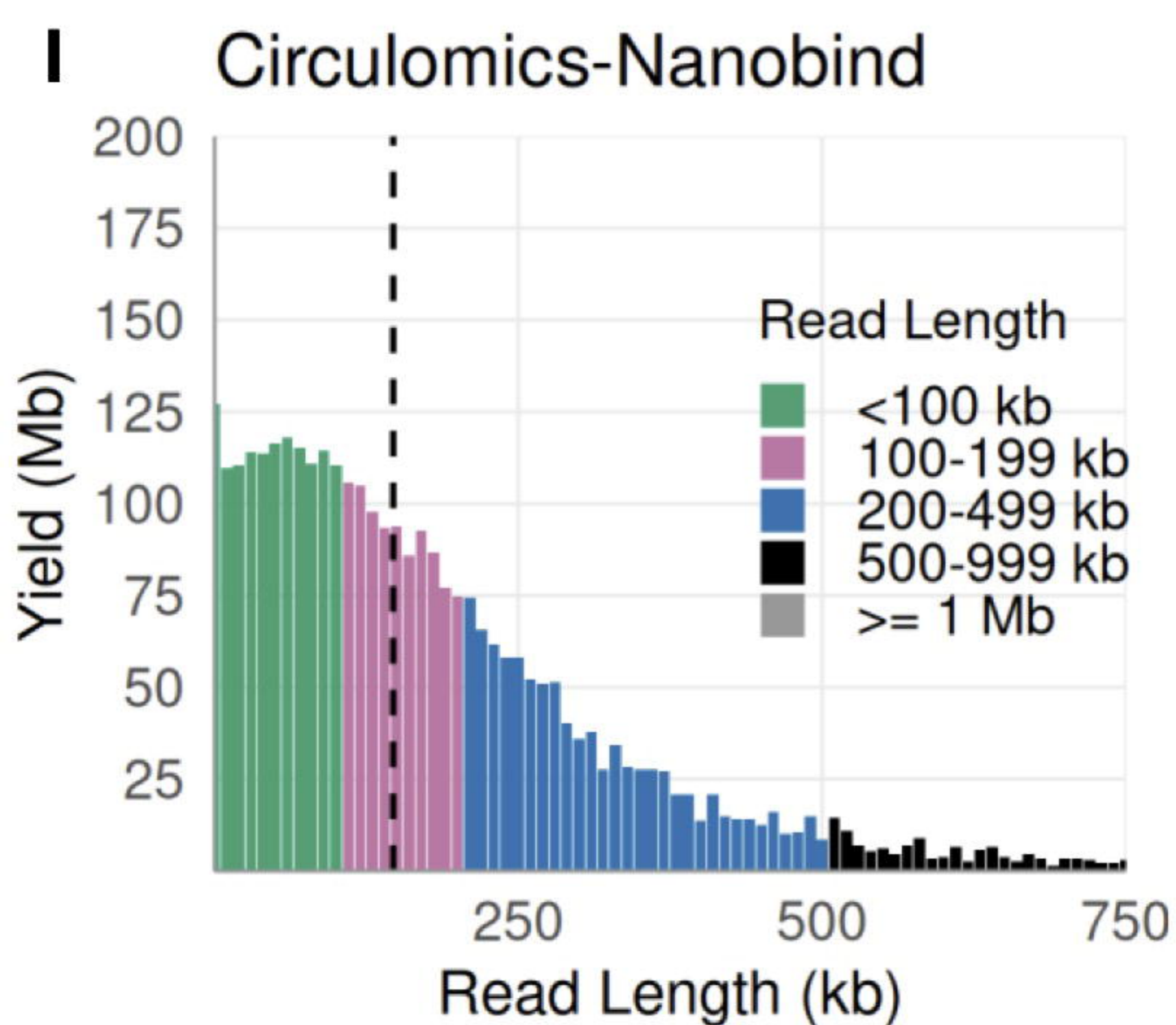

J


	Extraction-Clean up	N50 (bp)	Pores Used	Yield (% Bases $\geq$ 100 kb)	Average Occupancy (%)	Yield per Pore (Mb)
A	Phenol-Nemo	131,338	643	2.35 Gb (62.15%)	95.65	3.66
B	Phenol-Nemo 2	129,521	437	0.67 Gb (62.00%)	76.69	1.53
C	CTAB Lysis-Nemo	74,909	581	2.13 Gb (34.90%)	93.60	3.67
D	CTAB Lysis-Nemo 2	122,223	537	1.01 Gb (58.75%)	78.09	1.87
E	SDS Lysis-Nemo	106,027	461	0.92 Gb (52.55%)	76.48	1.99
F	CTAB and CoHex Lysis-Nemo	83,819	834	2.41 Gb (39.62%)	77.30	2.89
G	Monarch 300 rpm-Nemo	124,929	775	2.10 Gb (60.18%)	94.43	2.71
H	Monarch 900 rpm-Nanobind	106,493	862	3.35 Gb (52.70%)	97.36	3.89
I	Circulomics-Nanobind	147,123	1025	3.27 Gb (64.85%)	95.01	3.19

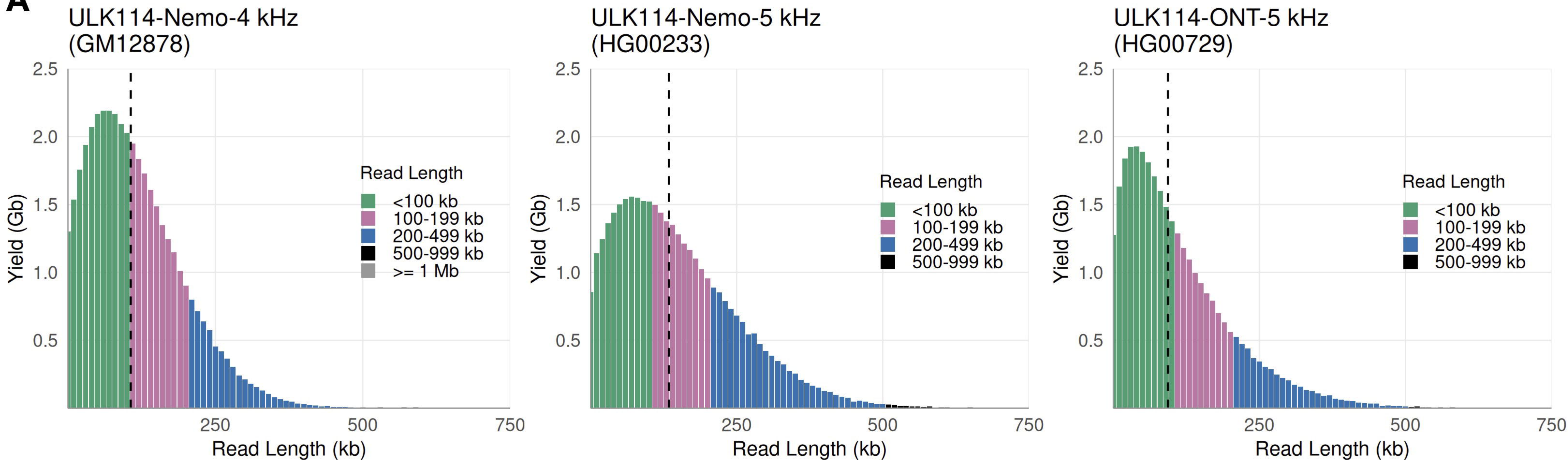A**B**


Sequencing metrics at 24 hours of run on PromethION platform

Sample	N50 (bp)	Pores Used	Yield (Reads >100 kb)	Yield per Pore (Mb)	Average Occupancy (%)	Mapping Frequency (Primary / Total Mapped)
ULK114-Nemo-4 kHz	102,692	7817	38.04 Gb (51.26%)	4.87	98.70	92.1% / 94.8%
ULK114-Nemo-5 kHz	133,881	7867	37.43 Gb (63.36%)	4.76	98.65	96.9% / 98.3%
ULK114-ONT-5 kHz	93,538	7879	32.30 Gb (47.08%)	4.10	98.74	96.9% / 98.0 %

A – PromethION Platform

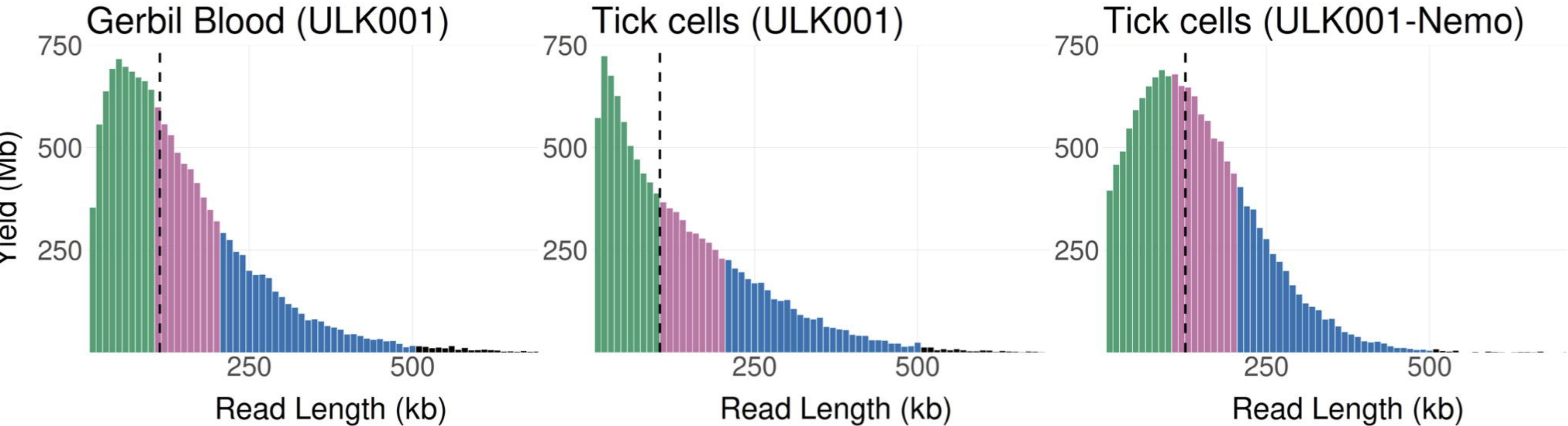


B – MinION Platform

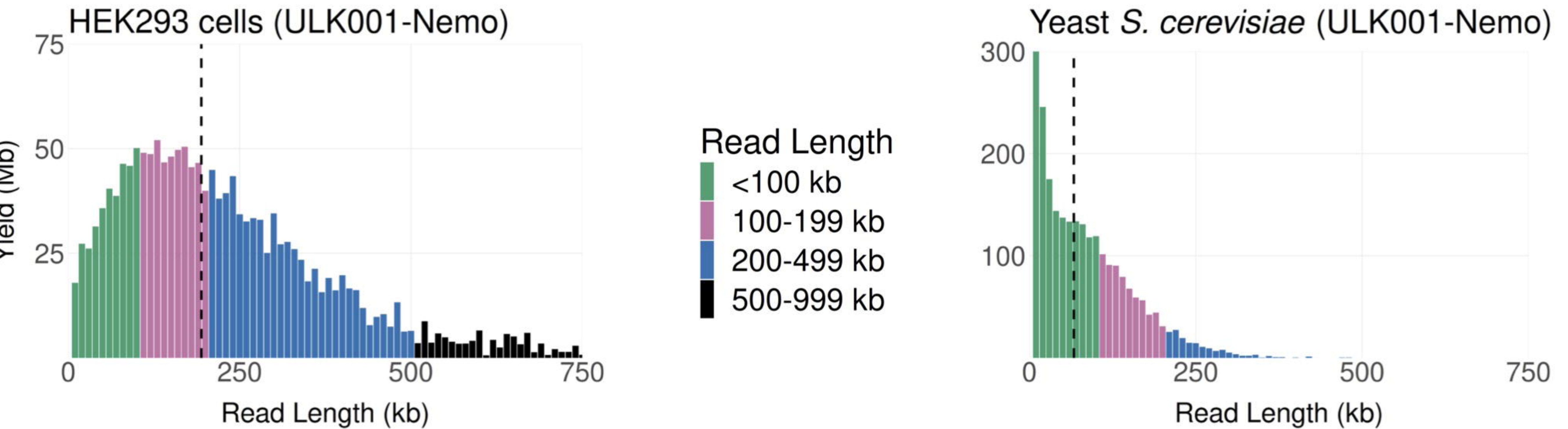

C – Sequencing Metrics at 6 Hours of Run

Sample (Library Protocol)	Platform	N50 (bp)	Pores Used	Yield (% Bases ≥ 100 kb)	Yield per Pore (Mb)	Average Occupancy (%)
Gerbil blood cells (ULK001)	PromethION	112,940	6220	14.16 Gb (55.43%)	2.28	79.89
Tick cells (ULK001)	PromethION	104,993	3492	11.11 Gb (51.64%)	3.18	87.50
Tick cells (ULK001-Nemo)	PromethION	126,008	3902	15.01 Gb (61.46%)	3.85	88.20
HEK293 cells (ULK001-Nemo)	MinION	194,200	728	1.63 Gb (77.94%)	2.24	91.53
Yeast <i>S. cerevisiae</i> (ULK001-Nemo)	MinION	66,694	762	2.47 Gb (33.36%)	3.24	92.87

A

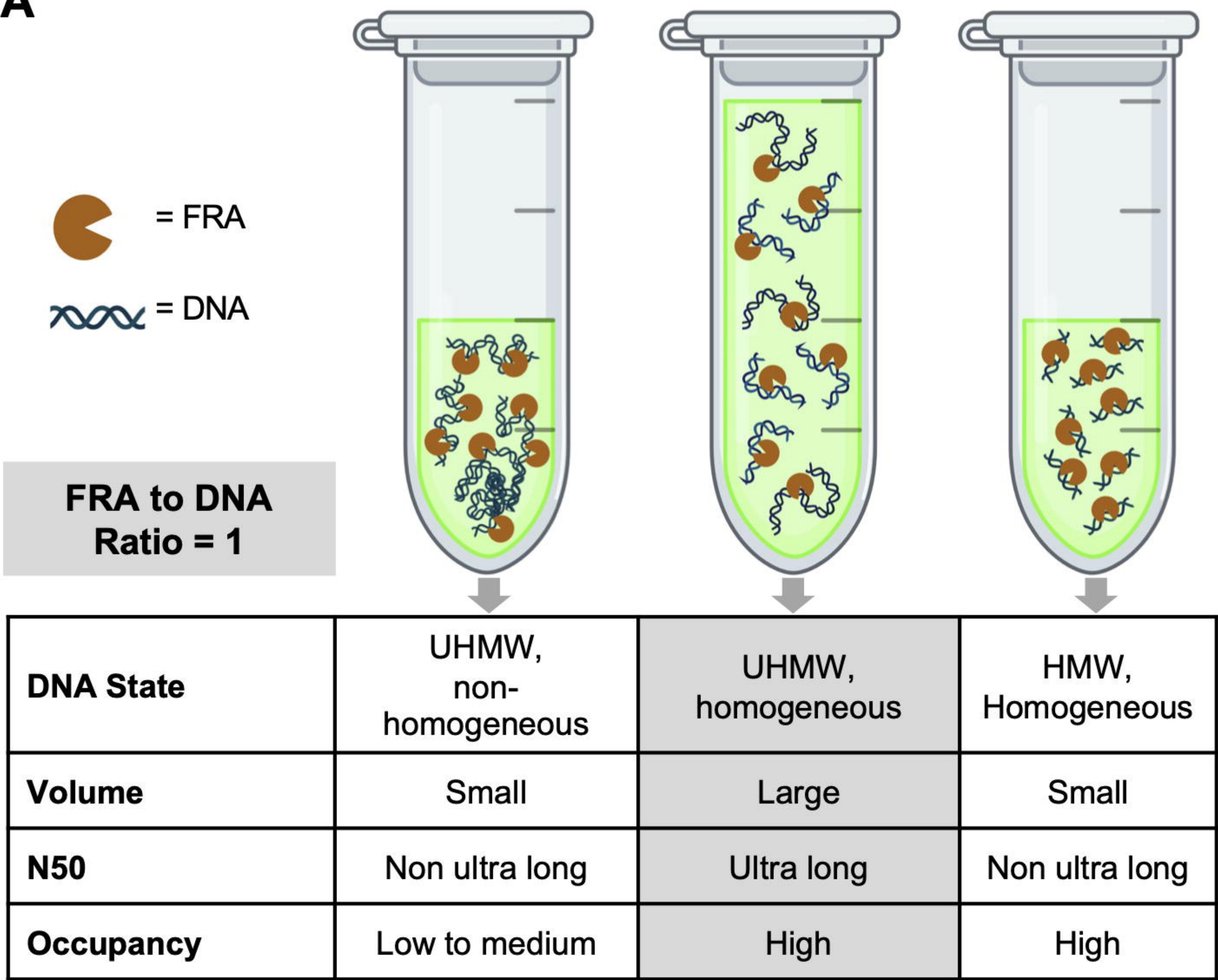

B

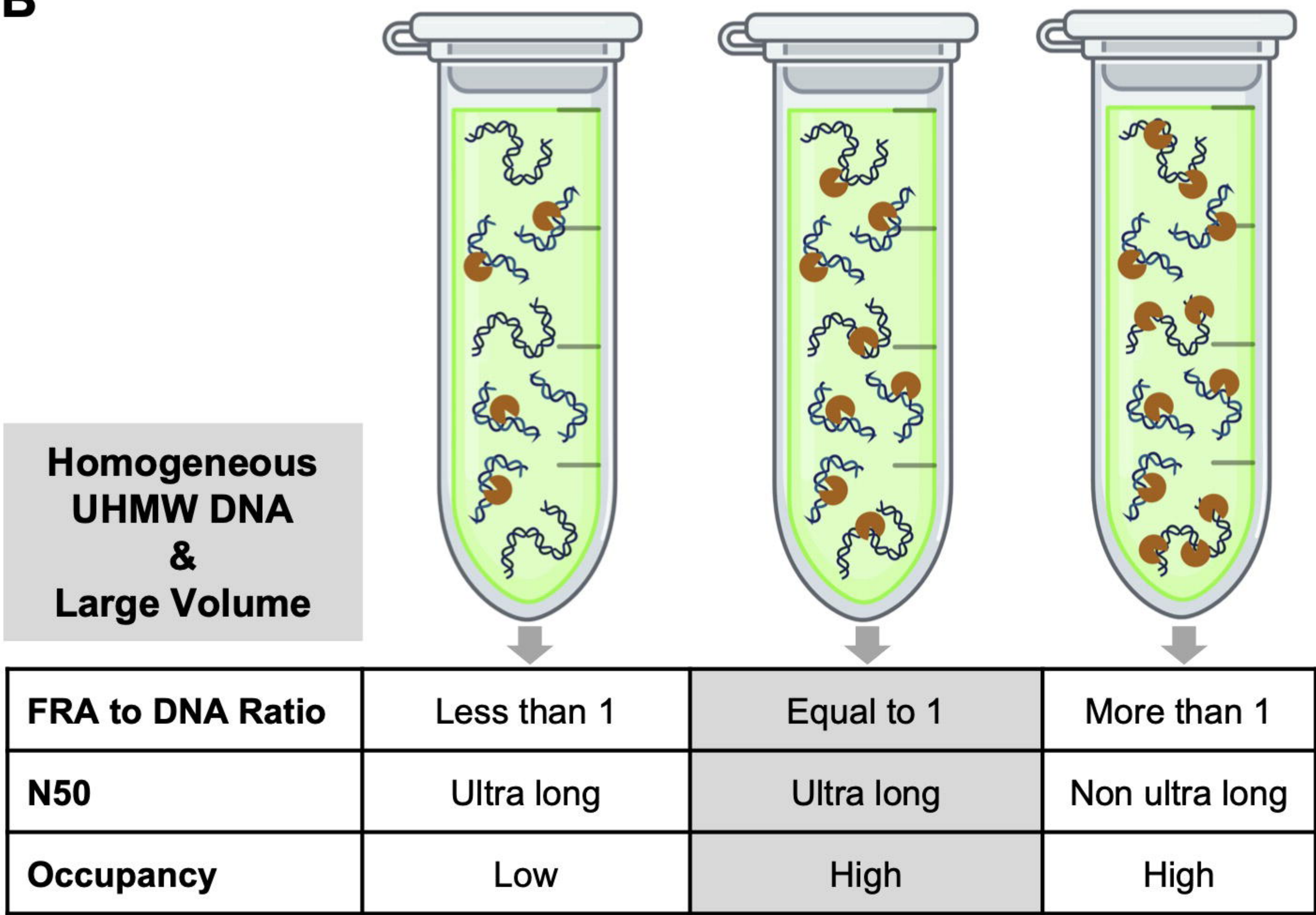


Table 1 Compilation of ultra-long (UL) library preparation protocols, optimized in human GM12878 cells

	Previous Protocols				New Protocols				
	Quick's UL (Quick, 2018)	Logsdon's UL (Logsdon, 2020)	Rocky Mountain (Tyson, 2019)	Krazy StarFish (KSF)	ULK001	RAD004- Nemo	ULK001- Nemo	ULK114	ULK114- Nemo
	A	B	C	D	E	F	G	H	I
Input cell number*	50 million	20-70 million	>2 million	>2 million	6 million	1-6 million	1-6 million	6 million	6 million
Input DNA amount	15 µg	2-3 µg	≥7.5 µg	≥7.5 µg	40 µg	5-40 µg	5-40 µg	40 µg	40 µg
ONT Library preparation kit	RAD004	RAD004	LSK109	RAD004	ULK001	RAD004	ULK001	ULK114	ULK114
Fragmentation Volume	Small (15 µl)	Small (18 µl)	Medium (60-100 µl)	Medium (2 × 100 µl)	Large (1 ml)	Large (0.5-1 ml)	Large (0.5-1 ml)	Large (1 ml)	Large (1-1.2 ml)
Library clean-up protocol	No	No	No	Partial using filter paper	Nanobind (Circulomics)	Nemo	Nemo	No	Nemo
Processing time (excluding incubations)	15 min	2.5 hours	18-24 hours	30 min	~2.5 hours	~2.5 hours	~2.5 hours	~2.5 hours	~2.5 hours
Input library per MinION flow cell	10-15 µg	2-3 µg	1 to 3 µg	2 to 4 µg	5 to 6 µg	1 to 5 µg	1 to 5 µg	NA	NA
No. of loads per library (MinION)	1	1	1 to 3	3	6	1 to 6	1 to 6	NA	NA
Expected yield per MinION flow cell	~3 Gb	1-2 Gb	10 to >20 Gb (depends on N50)	6-8 Gb (3 loads)	>20 Gb (3 loads)	>15 Gb (3 loads)	~20 Gb (3 loads)	NA	NA
Input library per PromethION flow cell	NA	NA	NA	NA	10-13 µg	10-13 µg	10-13 µg	10-13 µg	10-13 µg
No. of loads per library (PromethION)	NA	NA	NA	NA	3	1 to 3	1 to 3	3	3 to 4
Expected yield per PromethION flow cell (Gb)	NA	NA	NA	NA	NA	NA	24 Gb (1 load)	~85 Gb	up to >100 Gb
Ultra-long read length (≥100 kb)	Yes	Yes	No	Yes	Yes	Yes	Yes	Yes	Yes

(*) cell number was dependent on the extraction protocol used; here it is listed as the number of cells required by the original protocol; if cell number is unknown, use the DNA amount; (NA) not available

Table 2 Compilation of UHMW DNA extraction methods for UL sequencing, optimised in human GM12878 cells

	Without Kits (Phenol-based)			Without Kits (Phenol-free)		Kit-based	
	Quick's Phenol Protocol (Quick, 2018)	Logsdon's Phenol Protocol (Logsdon, 2021)	Phenol + Glass beads	CTAB or SDS Lysis + glass beads	CTAB-CoHex Lysis + glass beads	NEB Monarch Cell & Blood Kit	Circulomics Nanobind CBB Big DNA Kit
	A	B	C	D	E	F	G
Standard cell input	50 million	20-70 million	1-5 million	1-5 million	1-5 million	0.5-5 million	6 million
Phenol-chloroform	Yes	Yes	Yes	None	None	None	None
Lysis buffer	SDS	SDS	SDS	CTAB or SDS	CTAB and CoHex	Proprietary	Proprietary
SPRI Substrate	None	None	3-mm glass beads	3-mm glass beads	3-mm glass beads	4-mm glass beads	Nanobind disc
Precipitation reagents	Ethanol, Ammonium Acetate	Ethanol, Ammonium Acetate	Ethanol, Ammonium Acetate	Isopropanol, Ammonium Acetate	CoHex	Isopropanol, Proprietary binding buffer	Isopropanol, Proprietary binding buffer
DNA extraction time	4 hours	4 hours + overnight	2 hours	1 hour	1 hour	~30 min	2 hours
DNA resuspension time	2 days	2 days	min. overnight	min. 2 hours	min. overnight	min. 1 hour (standard overnight)	overnight

Table 3. Ratio of FRA to DNA in different genomes

Species	Genome Size	Cell Input Number	DNA Equivalent (μg)	FRA Volume (μl)
Human GM12878	3.2 Gb	1×10^6	5-6	1
Human HEK293	3.2 Gb	1×10^6	6-8*	1
Mongolian Gerbil	~2.7 Gb	Unknown	10	1.5
Tick <i>A. variegatum</i>	~6.0 Gb	1×10^6	11-13	2.5
Yeast <i>S. cerevisiae</i> (n)	12 Mb	2×10^8	4-5	6

(*) extracted DNA yield is likely to be more than normal diploid human cell line because of its hypotriploid nature (Lin 2014)