




Cell Load

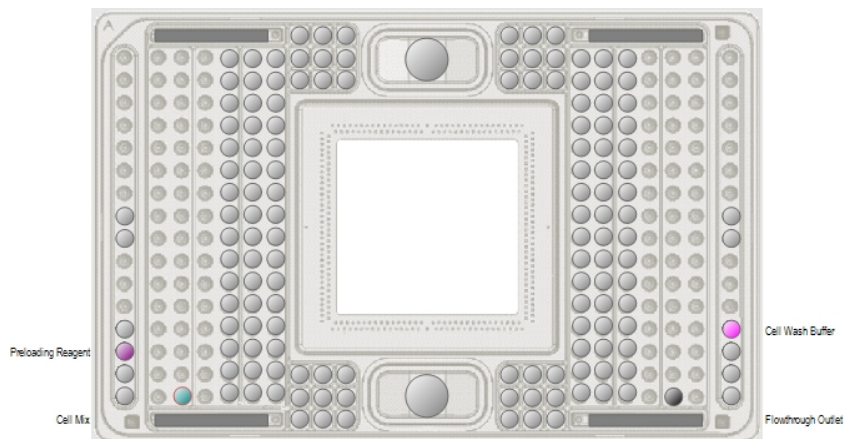

Script Run Time

Barcode	Estimate
1861x (5-10 um diameter cells)	0 hours, 15 minutes
1862x (10-17 um diameter cells)	0 hours, 30 minutes
1863x (17-25 um diameter cells)	0 hours, 27 minutes

Inlet Reuse

Name	Chip Inlet	Instructions
● Cell Mix	C1	Aspirate inlet prior to loading reagents
● Flowthrough Outlet	C2	Aspirate inlet prior to loading reagents (1862x, 1863x only)

Reagent Loading

Name	Volume	Chip Inlet	Notes
● Preloading Reagent	20 µl	2	
● Cell Wash Buffer	7 µl	5	
● Cell Mix	5 µl	C1	